Onduo Mobile Application

User Guide

Welcome to the Onduo mobile app, part of the **Onduo Virtual Care Software Platform!** Our mission is to support you on your journey to better manage Type 1 diabetes, Type 2 diabetes, hypertension, obstructive sleep apnea, and your mental and behavioral health. This smartphone application may help you:

- Track your glucose readings, blood pressure, and notice trends (Type 2 diabetes only)
 - Log your meals, exercise, and medications
 - Learn more about how food, exercise and other lifestyle activities impact your readings
 - Build and maintain healthy sleep behaviors and routines
 - Adjust to positive airway pressure (PAP) therapy and support ongoing therapy (Sleep Apnea users only)
 - Track your moods, mental health, and symptoms of depression and anxiety
 - Use self-help tools to make healthy changes and improve your mental well-being
 - Message directly with your Care Team, and access therapy or psychiatric care as needed (Mental & Behavioral Health users only)
 - Build healthy habits and skills
 - Connect with your Care Team for support
-

1.0 Directions for Use	10
1.1 Download and sign in to the app	12
1.2 Introduction to the app	13
1.3 Home tab: Your journey and activities	14
1.4 History tab: Exploring the impact on your glucose, blood pressure, PAP therapy, and mental and behavioral health	14
1.4.1 Symbols in your history	16
1.4.2 History with a blood glucose meter (BGM) - Type 2 diabetes only	17
1.4.3 History with a continuous glucose monitor (CGM) - Type 2 diabetes only	18
1.4.4 History with a blood pressure monitor	19

1.4.5 History with a Positive Airway Pressure (PAP) device	20
1.5 Connect your meter to the Onduo app	21
2.0 Device Description	22
3.0 Intended Use, Warnings, and Precautions	22
3.1 Intended use	22
3.2 Indications for use	22
3.3 Contraindications	23
3.4 Target population	23
3.5 Warnings	23
3.6 Precautions	23
3.7 Connected devices	25
3.8 Call, message, data, and cellular rates	26
3.9 App accessibility	27
3.10 Requirements	28
4.0 Your Care Team	28
4.1 Message and call your Care Team	29
4.1.1 Message your Care Team	29
4.1.2 Call your Care Team	30
5.0 Blood glucose meter (BGM) and manual readings (Type 2 diabetes only)	31
5.1. Intended use/Target population	31
5.1.1 Intended use	31
5.1.2 Target population	32
5.2 What's a BGM?	32
5.3 Connect your BGM	32
5.4 Manually log a blood glucose reading	33
5.5 Adjust your target glucose range	33
5.6 Review your blood glucose readings	34
5.7 Manage your blood glucose meter (BGM)	35
6.0 Continuous glucose monitor (CGM) (Type 2 diabetes only)	37
6.1 What's a CGM?	37
6.2 Connect your Dexcom CGM	37
6.2.1 Set up and insert your Dexcom G5 or G6 CGM	38
6.2.2 Connect your CGM data to the Onduo app	38
6.3 Track and adjust your daily goal	38
6.3.1 Adjust your time in range daily goal and your target range	39
6.4 Review your glucose history	39
6.5 Manage your Dexcom continuous glucose monitor (CGM)	40
7.0 Home sleep test (HST)	41

8.0 PAP device	42
8.1 What's a PAP device?	42
8.2 Set up and connect your PAP device	42
8.3 Using the PAP device	42
8.4 Review your PAP use	43
9.0 Medications	44
9.1 Add your medications	44
9.2 Edit and delete medications	45
9.3 Importing medications and medical supplies prescribed to you	45
9.4 Confirm medications imported from Surescripts or added by your Care Team	45
9.5 Log medications you have taken	46
10.0 Food, exercise, weight, and mood logging	47
10.1 Log food or a meal	47
10.2 Log exercise	47
10.3 Log mood	48
10.4 Log weight	48
10.5 View your logged events	49
10.6 Food logging suggestion	49
11.0 Step tracking	50
11.1 Track steps and adjust your step goal	50
11.2 View your steps	50
12.0 Your Journey: Activity sets, activities, checkpoints, insights, and education	52
12.1 Activity Sets	52
12.2 Learn about activities	52
12.3 Learn about insights and educational materials	52
12.4 Reminders for activities	53
12.5 Checkpoints	53
13.0 Video visits	54
13.1 Schedule a visit	54
13.2 (Old flow - No longer in use) Schedule a visit	54
13.3 (Old flow - No longer in use) View and modify your appointments	54
13.4 (Old flow - No longer in use) Join a video visit	55
14.0 Blood pressure	56
14.1 Indications for use/Target population	56
14.2 What's a blood pressure monitor?	57
14.3 Connect your blood pressure monitor	57
14.4 Manually log a blood pressure reading	57
14.5 Review your blood pressure and pulse readings	58

14.6 Finding your baseline blood pressure	59
14.7 Manage your blood pressure monitor	60
14.8 Blood pressure and pulse reading alerts	60
15.0 Obstructive Sleep Apnea	61
15.1 Indications for Use/Target population	61
15.2 Is Onduo's Sleep Apnea program right for me?	62
15.3 Getting an OSA diagnosis	62
15.4 Welcome Kit and PAP Therapy set up	65
15.5 Using PAP therapy	67
15.5.1 Schedule first night and first week activities	68
15.5.2 PAP therapy tracking	68
15.6 Using the Onduo app	69
15.7 Your guided learning activities	70
15.8 Messaging with your Care Team	70
15.9 Ordering replacement parts for your PAP device	70
16.0 Mental and Behavioral Health	71
16.1 Indications for use/Target population	71
16.2 Is Onduo's Mental and Behavioral Health right for me?	71
16.3 Getting started	72
16.4 Your journey	76
16.5 Messaging with your Care Team	76
16.6 Mental Health Checkpoints	77
16.7 Video visits with clinicians	78
17.0 Type 1 Diabetes	79
17.1 Intended use	79
17.2 Target population	79
17.3 Getting started	80
17.4 Devices and supplies	80
17.5 Video visits with Onduo PC	80
17.6 Diabetes device data dashboard	80
18.0 Settings	82
18.1 Set and change your schedule	82
18.2 Add and confirm medications	82
18.3 Adjust your step goal	82
18.4 (Old flow - No longer in use) My Appointments	83
18.5 Add and view connected meters	83
18.6 Set up and manage notifications	83
18.7 Change your target range and time in range (Type 2 diabetes only)	84

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
1.0 Directions for Use

PRECAUTIONS

IN THE EVENT OF A MEDICAL EMERGENCY, CALL 911 OR GO TO THE NEAREST EMERGENCY ROOM. The Onduo app is not intended to be used to respond to a medical emergency.

IF YOU ARE HAVING A MENTAL HEALTH EMERGENCY, CALL 911, GO TO THE NEAREST EMERGENCY ROOM, OR CALL/TEXT A CRISIS HOTLINE. If you are thinking about hurting yourself or others, or if you need immediate support, call 911 or have a friend or family member take you to the nearest emergency room. For confidential mental health support, you can also call the National Suicide Prevention Lifeline 24/7 at 1-800-273-8255 or text the Crisis Text Line (741-741). If you need immediate help, call 911 or have a friend or family member take you to the nearest emergency room. This guidance is also available for reference at the bottom of the Home tab:

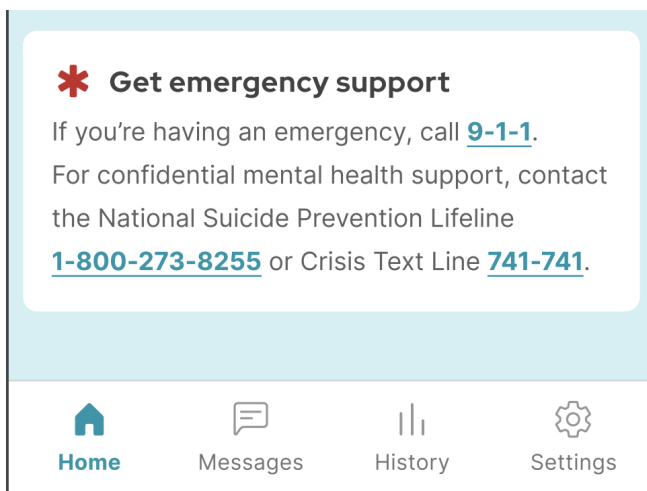Follow your physician's or mental health provider's advice if it conflicts with the app. The Onduo app is not intended as a substitute for your physician's guidance. Before making major lifestyle changes, please consult your physician. You should not make any treatment decisions or any medication adjustments based on data from the Onduo app.

This system is not intended to be used as a medication intake reminder system. Always defer to your prescribing physician's instructions. If you receive a medication logging notification at the wrong time, edit your medications list (refer to **Section 6.2**) and notify your Care Team.

Defer to your glucose device for primary glucose readings. If you believe the Onduo app is displaying incorrect glucose data, please defer to the readings on your glucose device. Please contact the Onduo support team if you believe the app is displaying incorrect data. Please note that only the Type 2 diabetes program supports in-app glucose data.

Defer to your blood pressure device for primary blood pressure readings. If you believe the Onduo app is displaying incorrect blood pressure data, please defer to the readings on your blood pressure device. Please contact the Onduo support team if you believe the app is displaying incorrect data.

Defer to your PAP device for primary PAP data. Follow up with your sleep physician if you have symptoms of concern related to your OSA or PAP therapy. Daytime sleepiness can continue even with regular PAP use. Never drive, operate heavy machinery, or perform other potentially risky tasks while experiencing sleepiness. The OSA Program includes scheduled visits with a sleep physician, and you may contact the physician outside of these visits if needed.

If you take insulin, consult with your physician before making any lifestyle changes. You may need to modify the amount of insulin you take if you change your lifestyle, including food and exercise.

Exercise can lower blood glucose levels, and in combination with certain medications, exercise may lead to hypoglycemia.

Versions:

iOS 14 or newer

Android version 8 or higher

Google Chrome browser version 96.0.4664.55 or higher

1.1 Download and sign in to the app

This IFU is for the latest version of the Onduo app. Please upgrade to the latest version if you haven't already, before reading and using this IFU. You can download the Onduo app in the [App Store](#) or [Google Play Store](#). Please note that if you are using a version older than the minimum supported version, you will be required to update your app in order to continue using it. When you open the app, you will see a link to the app store, or instructions on how to update your app.

You will need to log in to the Onduo app with the same account you used to register for the Onduo Program. The app allows you to sign in via the following two methods:

Google Account: If you used the Sign in with Google option, use this to sign in to the Onduo app. This account is only required to sign in to the app. The app does not associate any information with your Google Account.

Onduo Account: If you used your email to create a new Onduo account, use this to sign in to the Onduo app.

After logging in, if configured, you will see a page with Onduo and sponsor logos. Then you'll begin your experience in the app.

1.2 Introduction to the app

The app has 4 tabs:

-  **Home:** Choose your journey, work on activities, track your progress, and view educational content.
-  **Messages:** Message with your Care Team.
-  **History:** Track your data around glucose and blood pressure readings (Type 2 diabetes only), food/drink logs, movement logs (steps and other exercise), medication logs, mood logs, and PAP therapy data and insights (only available for users enrolled in Sleep Apnea). Review your depression score (only available for users enrolled in Mental and Behavioral Health).
-  **Settings:** Adjust the schedule of your daily routine, add/remove medications, adjust your daily goals, track your steps, edit your emergency contact information, and adjust other account settings.

Tap the task menu  button on the upper right corner of the Home tab to:

- Log a meal
- Track your exercise
- Log taking a medication
- Log a blood glucose reading (only available for users enrolled in Type 2 diabetes management)
- Log your mood
- Log a blood pressure reading (only available for users enrolled in hypertension management)

1.3 Home tab: Your journey and activities

The Home tab is where you can choose a health journey that's important to you, or see activity sets assigned to you. Your journey will include several activity sets for you to work on. Each activity set contains up to eight activities. Activities can include new skills to build, information to learn, or challenges to conquer. As you complete each activity, you'll unlock a new one. You'll also have the option to skip activities that don't work for you.

If you see journeys on the Home tab, you can click "View full journey" to view your progress and see what's coming up next. When you've completed a full journey, you'll be able to choose a new one. You may also see recommended educational content and activities in your Home tab, along with surveys from time to time. These surveys help us to better understand your health and personalize your program.


1.4 History tab: Exploring the impact on your glucose, blood pressure, PAP therapy, and mental and behavioral health

Your History tab displays information that may help you evaluate how your lifestyle impacts your glucose, blood pressure readings, and mental well-being. Your History tab displays information including:

- Glucose readings (in-app glucose readings only available for users enrolled Type 2 diabetes management)
- Blood pressure readings (only available for users enrolled in hypertension management)
- Food/drink
- Movement, including steps and other exercise
- Medication
- Mood
- PAP usage, AHI, and Mask Leak (only available for users enrolled in Sleep Apnea)
- Depression score (only available for users enrolled in Mental and Behavioral Health)

There are three ways to review your data:

- **Insights:** View data from how things like logged meals and exercise impact your blood glucose and blood pressure
- **Logbook:** View data categorized by blood glucose (Type 2 diabetes only), blood pressure, food/drink, movement, medication, mood, etc. Each subsection is broken into detailed views.
 - **Day:** View all your data logged on each day. Tap the arrows < > to see data for different days.
 - **Week:** View your data over the course of the week. Tap the arrows < > to see data for different weeks.
 - **Month:** View your data over the course of the month. Tap the arrows < > to see data for different months (you'll only see daily and weekly data if you're an OSA member).
 - **Checkpoints:** View completed checkpoints as well as when you started taking checkpoints with us.

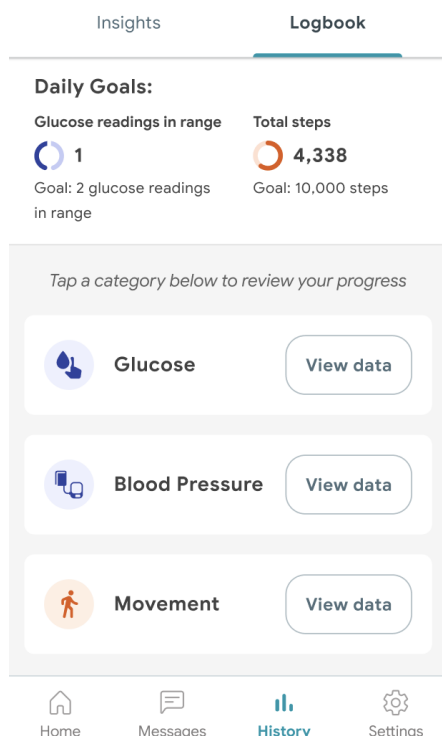1.4.1 Symbols in your history


See the table below to understand what the symbols in your history mean.

Symbol	What it means
	Daily goal not met (Type 2 diabetes only): This symbol is shown in your day and week history when you have not yet met (or did not meet) your daily target range or PAP usage goal for that day.
	Daily goal met (Type 2 diabetes only): This symbol is shown in your day and week history and Home tab when you have met your daily target glucose range or PAP usage goal for that day.
	Target glucose range (Type 2 diabetes only): This blue bar shows your target glucose range. If a reading is inside this blue bar, that means a reading is in-range. If it is above or below the blue bar, that means that a reading is out of range.
	Reading in-range (Type 2 diabetes only): This symbol shows blood glucose readings that you have taken that are within your target range.
	Reading out of range (Type 2 diabetes only): This symbol shows blood glucose readings that you have taken that are not within your target range.
	Meal: This symbol shows when you logged a meal.
	Exercise: This symbol shows when you logged exercise.
	Medication: This symbol shows when you logged that you took medication.
	Mood: This symbol shows when you have logged your mood.

	Steps: These vertical lines show the steps you took over the course of the day.
	Survey: This symbol shows your survey results.
	AHI: This symbol shows your AHI, or apnea-hypopnea index. AHI measures the average number of abnormal breathing events detected per hour of PAP use.
	Mask leak: This symbol shows your mask leak level. Mask leak is an estimate of the total rate of air unintentionally escaping your mouth or mask.
	Average PAP usage: This symbol shows your average amount of time you've used your PAP device over 2 or more nights.

1.4.2 History with a blood glucose meter (BGM) - Type 2 diabetes only

For those using a **BGM** to manage type 2 diabetes, your Glucose category of the History tab will show your logged events displayed in the context of your glucose readings. Please note that in-app glucose data is only supported in the Type 2 diabetes management program.

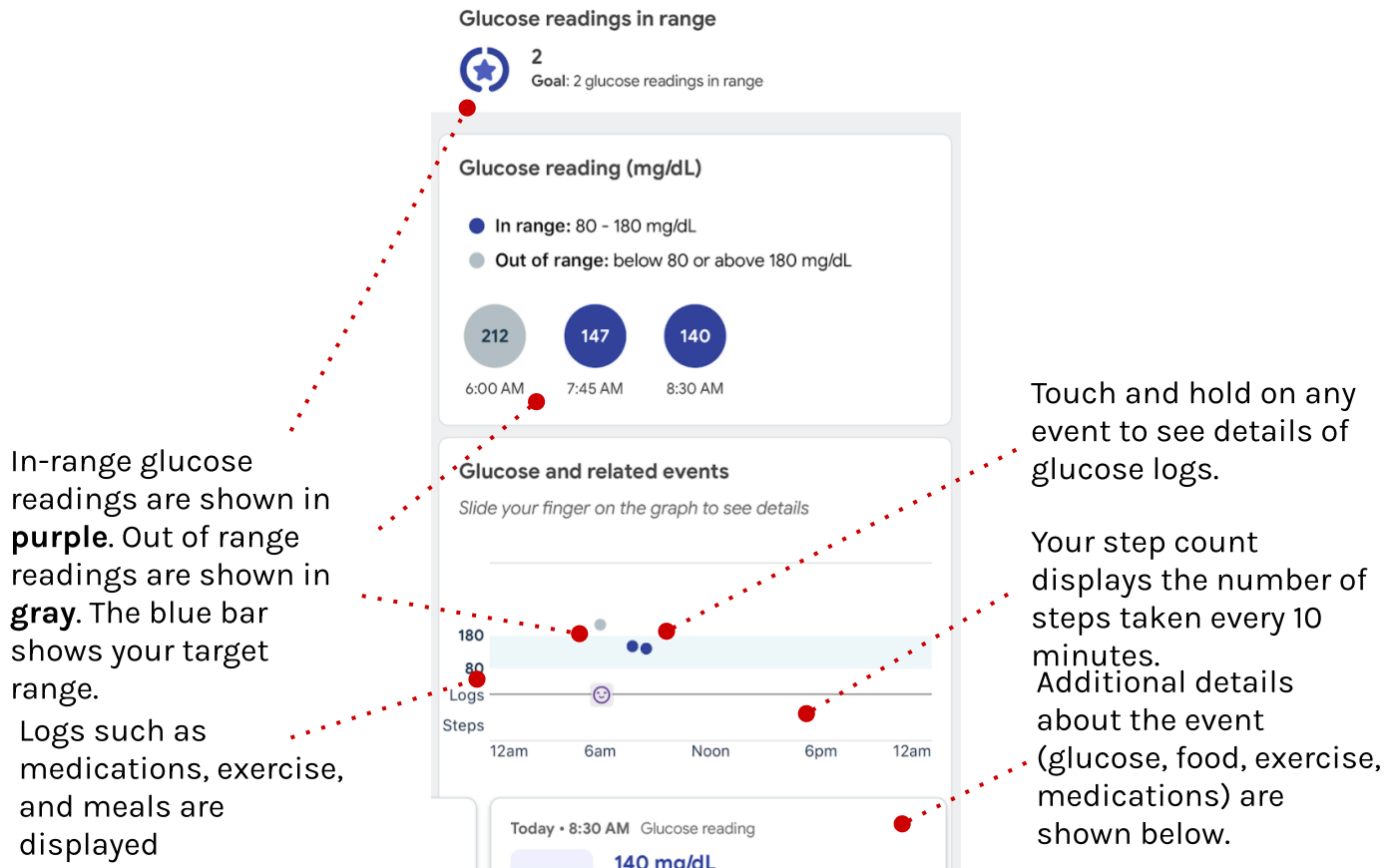Scroll left and right through the cards on the bottom of the page to see information about each day. Tap a selected card to show more information.

1.4.3 History with a continuous glucose monitor (CGM) - Type 2 diabetes only

For those using a **CGM** to manage type 2 diabetes, your Glucose category of the History tab will show your logged events displayed on your glucose graph.

See the progress you've made on your time in range goal for the day. In-range glucose readings are shown in **purple**. Out-of-range readings are shown in **gray**. Readings from your blood glucose meter will be shown as dots on your graph.

Touch and hold on any event to see details of glucose data.

Your step count displays the number of steps taken every 10 minutes.

Additional details are shown below the graph or in these cards.

Logs such as medications, exercise, and meals are displayed

Scroll left and right through the cards on the bottom of the page to see information about each day. Tap a selected card to show more information.

1.4.4 History with a blood pressure monitor

For those enrolled in the hypertension management program, your Blood Pressure category of the History tab will show blood pressure and pulse readings from your connected blood pressure monitor. Manually logged readings will also show up in this table.


Blood Pressure

[What do blood pressure and pulse numbers mean?](#)

Date	Time	SYS/DIA (mmHg)	Pulse (BPM)
Tue, Feb 11	10:24 PM	81/54	65
Tue, Feb 11	7:20 AM	116/72	143
Tue, Feb 11	6:33 AM	130/105	76

1.4.5 History with a Positive Airway Pressure (PAP) device

For those enrolled in the obstructive sleep apnea program, the PAP Therapy category in your History tab will show PAP usage, Apnea Hypopnea Index (AHI) and Mask Leak data from your connected PAP device.

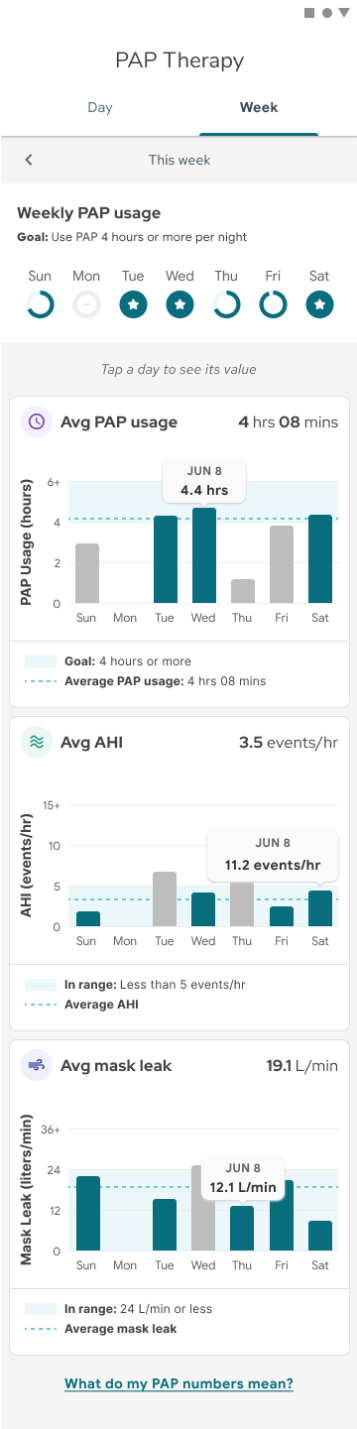1.5 Connect your meter to the Onduo app

1. Open the Onduo app.
2. Follow the instructions in the app to connect your app with your connected meter.

Helpful tips:

- For instructions, open the Onduo app and look on the Home tab or under **Settings > Connected meters**.

2.0 Device Description

The Onduo Software Platform acquires health data through various sources based on the users health condition such as commercially available medical devices like FDA-cleared blood glucose meters or FDA-approved continuous glucose monitoring systems (For Type 2 diabetes), or third party data pooler (For Type 1 diabetes), FDA-cleared blood pressure monitors (hypertension), FDA-cleared PAP devices and masks (obstructive sleep apnea), or self provided by users. The Platform also connects to commercially available non-medical fitness apps to receive activity data. Lastly, the Platform connects to Verily developed non-medical device software like the Registration server and Integration server that aid in operational support and administrative functions like adding new users to the Platform and enabling transfer of device data along with other pilot logistics data, respectively. Patient data can be provided in a Clinician Report.

3.0 Intended Use, Warnings, and Precautions

3.1 Intended use

The Onduo Software Platform is intended to provide individuals with health conditions, their healthcare professionals, wellness, and care team a review, analysis, and evaluation of health data using the platform tools and modules. Users are able to view the results for any analysis done by the software.

The Onduo software platform is not intended for real-time monitoring of health data, nor is it meant as a primary display of data. The Onduo Software Platform is not intended to provide treatment decisions nor is it to be used as a substitute for professional healthcare advice. The platform is meant to be used adjunct to a member's standard of care.

3.2 Indications for use

The Onduo Mobile App is intended to be used by adults aged 18 or older on iOS and Android mobile operating systems with one or more of the following conditions: diabetes, hypertension, obstructive sleep apnea, and/or mental or behavioral health concerns.

3.3 Contraindications

The Onduo Software Mobile App is not intended for use by individuals:

1. who have had an inpatient admission to the hospital for a mental health reason OR a suicide attempt, in the last 3 months
2. who are currently pregnant
3. who have ever been diagnosed with any of the following conditions: Severe or end-stage kidney disease, Cystic fibrosis, Cirrhosis or liver failure
4. who have had an organ transplant or bone marrow transplant

3.4 Target population

Health data surfaced through the Onduo platform may assist trained medical professionals in making healthcare decisions by trained medical professionals. The intended population for each disease condition module is described in detail for each of the disease conditions for each module as described in each of the sections below.

3.5 Warnings

The Onduo app is not intended as a substitute for your primary care physician's guidance. Before making major lifestyle changes, please discuss your treatment plan with your primary care provider.


You should not make any treatment decisions or any medication adjustments based on data from the Onduo app.

3.6 Precautions

IN THE EVENT OF A MEDICAL EMERGENCY, CALL 911 OR GO TO THE NEAREST EMERGENCY ROOM. The Onduo app is not intended to be used to respond to a medical emergency.

IN THE EVENT OF A MENTAL HEALTH EMERGENCY, CALL 911, GO TO THE NEAREST EMERGENCY ROOM, OR CALL/TEXT A CRISIS HOTLINE. If you are thinking about hurting yourself or others, or if you need immediate support, call 911 or have a friend or family member take you to the nearest

emergency room. For confidential mental health support, you can also call the National Suicide Prevention Lifeline 24/7 at 1-800-273-8255 or text the Crisis Text Line (741-741). This guidance is also available for reference at the bottom of the Home tab:

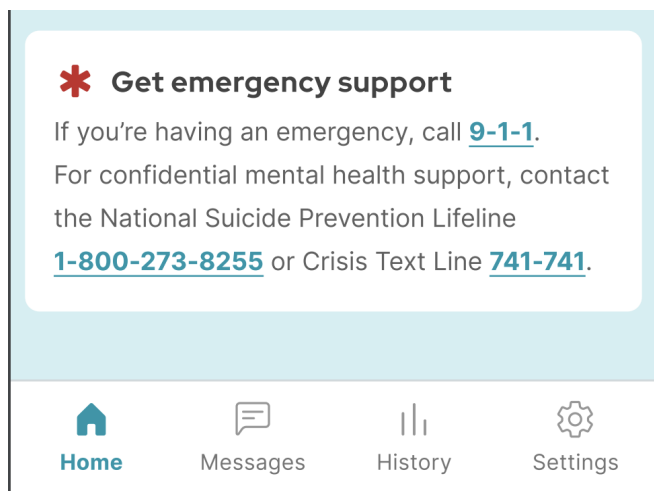Follow your physician's or mental health provider's advice if it conflicts with the app. The Onduo app is not intended as a substitute for your physician's guidance. Before making major lifestyle changes, please consult your physician.

If you attempt to use the Onduo app on an untested operating system, the experience may not work as intended or documented in this instructional guide. The Onduo app is designed to work with the following operating systems:

- The Android app is functional on Android OS 8.0 and higher
- The Apple app is functional on iOS 14.0 and higher

Defer to your glucose device for primary glucose readings. If you believe the Onduo app is displaying incorrect glucose data, please defer to the readings on your glucose device. Please contact the Onduo support team if you believe the app is displaying incorrect data. Please refer to the specific device's guide for complete warnings and safety information. Please note that only the Type 2 diabetes program supports in-app glucose data.

Glucose readings transferred from a CGM will be displayed on a 3 to 4-hour delay in the Onduo app.

Always consult with a physician before changing your target range. The American Diabetes Association (ADA) recommends that you maintain a glucose range of 80-180 mg/dL.

This system is not intended to be used as a medication intake reminder system. Always defer to your prescribing physician's instructions. If you receive a medication logging notification at the wrong time, edit your medications list (refer to **Section 6.2**) and notify your Care Team. Medication notifications remind you to log your medication intake. To turn on medication notifications refer to **Section 9.1** to adjust your notification settings

If you take insulin, consult with your physician before making any lifestyle changes. You may need to modify the amount of insulin you take if you change your lifestyle, including food and exercise.

Exercise can lower blood glucose levels, and in combination with certain medications, exercise may lead to hypoglycemia.

Sign out of the Onduo app to keep your data confidential. If you have not enabled a passcode on your phone, you will be automatically logged out of the Onduo app after ten minutes of inactivity in order to keep your data confidential. If you are logged out, you'll be directed to your phone's settings to set up a passcode.

3.7 Connected devices

The app can connect wirelessly to commercially available, FDA-regulated medical devices as well as non-medical devices through the participant's smartphone.

Depending on the program you enroll in, the following devices may be available or prescribed to you:

- FDA-cleared blood glucose meter (e.g.: Telcare BGM)
- FDA-approved CGM (e.g.: Dexcom G5® Mobile CGM System)
- FDA-cleared blood pressure monitor (e.g.: Omron HEM-9210T)
- FDA-cleared PAP machine and accessories
- FDA-cleared home sleep test (e.g.: WatchPAT™ ONE)
- A&D UA-651BLE Blood Pressure Monitor
- CVS Health Series 800 Blood Pressure Monitor
- Omron 7 Series (BP7350) Wireless Upper Arm Blood Pressure Monitor
- Omron 5 Series (BP7250) Upper Arm Blood Pressure Monitor
- Omron HEM-9200T Blood Pressure Monitor
- Omron HEM-9210T Blood Pressure Monitor
- [Omron BP9300T Blood Pressure Monitor](#)
- [Omron BP9310T Blood Pressure Monitor](#)
- Pyle Health PHBPB20 Blood Pressure Monitor
- Welch Allyn H-BP100SBP (Series 1700) Blood Pressure Monitor
- Welch Allyn RPM-BP100 (Series 1500) Blood Pressure Monitor
- Zewa UAM-820BT Blood Pressure Monitor

Please refer to the specific device's guide for complete warnings and safety information.

3.8 Call, message, data, and cellular rates

Regular usage of the app can happen over a Wi-Fi connection, but may use cellular data if no Wi-Fi is available.

Some examples of mobile data usage in the app may include:

- Uploading a photo with a meal
- Messaging your Care Team

Standard phone rates apply when calling your Care Team.

Overall, the Onduo app uses a small amount of data compared to the total used by most people in a day. The app will send this information via Wi-Fi if you are connected to a wireless network, otherwise it will use the data from your phone plan. You are responsible for all data usage charges for your plan.

Without cellular or Wi-Fi connectivity, you will not be able to send messages in the app and your Care Team will not be able to see newly uploaded glucose values.

There are a few things that will still work if you're offline (without access to a cellular or Wi-Fi connection):

- Logging blood glucose readings manually with a BGM (Type 2 diabetes only)
- Logging your medication
- Logging your mood
- Viewing your history
- The step counter
- Logging a meal

You won't be able to reach your Care Team or use other parts of the app. Information and messages will sync when your mobile device connects to a cellular or Wi-Fi network.

3.9 App accessibility

The Onduo app supports the use of assistive technologies to ensure users with disabilities have equal access to information and data. Assistive technologies include:

- **iOS:** VoiceOver, Zoom, Assistive Touch, Invert Colors, Speak Screen, Closed Captioning, Larger Text, Voice Access, and Switch Access.
- **Android:** TalkBack, Magnification, Closed Captioning, Large Font Size, Large Display Size, Voice Control, and Switch Control.

To enable assistive technology on your phone, go to Settings > Accessibility, and select the desired tool.

3.10 Requirements

The Onduo app is designed to work with the following requirements:

- Your phone must have a cellular or Wi-Fi connection
- For the best experience all Onduo app notifications should be enabled
- If you use an iPhone, you must be running iOS 14 or above
- If you use an iPhone, for the best experience the Onduo app should be granted access to HealthKit
- If you use an Android phone, you must be running Android 8.0 Oreo or above
- If you use an Android phone, for the best experience the Onduo app should be granted access to Google Fit
- If you use an Android phone, for the best experience, Bluetooth should be on and Bluetooth permissions should be granted to the Onduo app
- If you will be using a Dexcom G6 CGM, you must download the Dexcom G6 Mobile System app

4.0 Your Care Team

As part of using this app, you may be participating in a larger program with Onduo that gives you access to a Care Team.


PRECAUTION

IN THE EVENT OF A MEDICAL EMERGENCY, CALL 911 OR GO TO THE NEAREST EMERGENCY ROOM. The Onduo app is not intended to be used to respond to a medical emergency.

IF YOU ARE HAVING A MENTAL HEALTH EMERGENCY, CALL 911, GO TO THE NEAREST EMERGENCY ROOM, OR CALL/TEXT A CRISIS HOTLINE. If you are thinking about

hurting yourself or others, or if you need immediate support, call 911 or have a friend or family member take you to the nearest emergency room. For confidential mental health support, you can also call the National Suicide Prevention Lifeline 24/7 at 1-800-273-8255 or text the Crisis Text Line (741-741). This guidance is also available for reference at the bottom of the Home tab:

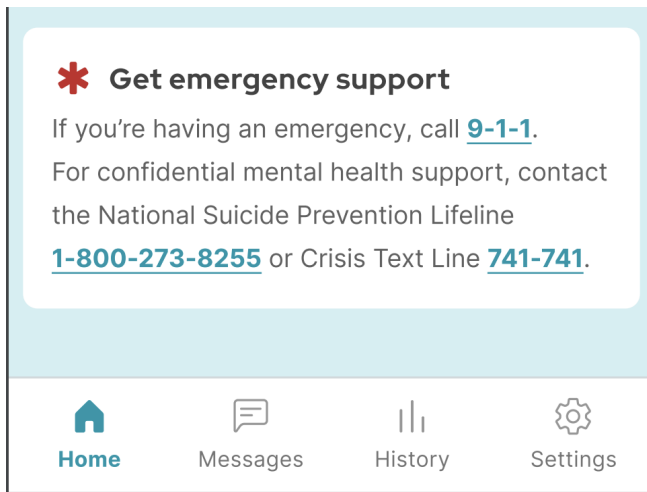Follow your physician's or mental health provider's advice if it conflicts with the app. The Onduo app is not intended as a substitute for your physician's guidance. Before making major lifestyle changes, please consult your physician.

4.1 Message and call your Care Team

4.1.1 Message your Care Team

You can message your Care Team through the Onduo app.

1. Open the Onduo app.
2. Tap **Messages** .
3. Tap the text box then start typing.
4. [Optional] If you want to take a picture to send to your Care Team, tap the Camera  icon to add a photo.

5. [Optional] If you want to send a picture you already took to your Care Team, tap Album .
6. When done, tap send  to message your Care Team.

Please note that member messaging and notes are not intended to be used as the official medical record for diagnosis and treatment. The communications for diagnosis and treatment must be communicated in the electronic medical record.

4.1.2 Call your Care Team

You may also be able to call your Care Team, though it's best to send a message to confirm that a Team member is available to talk.

1. Open the Onduo app.
2. Tap **Messages**  > Call , if you see this icon.

5.0 Blood glucose meter (BGM) and manual readings (Type 2 diabetes only)

Blood glucose meters (BGMs) are an important part of your experience with the Onduo app if you're using the Onduo app to help manage type 2 diabetes. We recommend you use the BGM provided to you. You can use your own meter if you prefer, but you may not be able to connect the device to the Onduo app and may have to log a reading manually. Please note that only the Type 2 diabetes program supports in-app glucose data.

Start by learning how to connect your meter and review your readings.

PRECAUTION

Defer to your glucose device for primary glucose readings. If you believe the Onduo app is displaying incorrect glucose data, please defer to the readings on your glucose device. Please contact the Onduo support team if you believe the app is displaying incorrect data. For specific instructions for using the device provided to you (such as strip insertion or on-device alerts), please refer to the meter's instructions for use.

5.1. Intended use/Target population

5.1.1 Intended use

The Onduo Diabetes Management Software Platform is intended to provide individuals with Type 2 diabetes and their healthcare professionals, wellness and care teams with a review, analysis and evaluation of retrospective glucose data from compatible FDA-cleared blood glucose meters and FDA-approved continuous glucose monitoring systems in support of an effective diabetes management program.

The Onduo Diabetes Management Software Platform is not intended to provide treatment decisions, nor is it to be used as a substitute for professional healthcare advice.

5.1.2 Target population

The Onduo Diabetes Management Software Mobile App is intended to be used by individuals who:

- Have been diagnosed with Type 2 Diabetes,
- Are aged 18 and older, and
- Are using iOS or Android mobile operating systems.

The Onduo Diabetes Management Care Team Console is intended to be used by Onduo healthcare professionals and care team users.

5.2 What's a BGM?

You may receive a blood glucose meter (BGM) as part of a larger Onduo program. A BGM is a device that measures the amount of glucose in your blood.

5.3 Connect your BGM

If you are part of a program in Onduo with a connected BGM, connect your meter to the Onduo app to automatically save and track your blood glucose readings in the Onduo app.

1. Open the Onduo app.
2. On the Home  tab, under "Connect your blood glucose meter," tap **CONNECT**.
If you don't see this, open the Onduo app **Settings > Connected meters > BioTel Care Meter**.
3. Follow the instructions in the app, including:
 - a. Locate the serial number on the back of your meter.
 - b. Verify the serial number of your meter in the Onduo app.
 - c. Program your name into the meter.

4. The app will confirm once your meter has been successfully connected.
5. The first time a blood glucose reading is received from your meter, you'll receive a "Celebration" push notification with a link to your Blood Glucose View Data screen. You'll also see Celebrations when your app receives blood glucose data for two or more days in a row, and when data is received after a break of seven days or more.

If you have questions, or if you don't see the option to "Connect your blood glucose meter," please contact your Care Team. Or, try opening the Onduo app **Settings > Connected meters > BioTel Care Meter**.

5.4 Manually log a blood glucose reading

You can log glucose readings manually:

1. Open the Onduo app.
2. Tap the task menu  **> Glucose**.
3. Type in your blood glucose reading (recorded as mg/dL).
4. Adjust the logged your glucose reading, if needed.
5. Tap the buttons to provide more information about the reading:
 - **Fasting:** Fasting means after not having anything to eat or drink (except water) for at least 8 hours before the reading.
 - **Bedtime:** This reading was taken at bedtime.
 - **Before meal:** Tap to specify if this reading was before a particular meal.
 - **After meal:** Tap to specify if this reading was after a particular meal.
6. Tap **Save**.
7. The first time you manually log a blood glucose reading, you'll see a "Celebration" message in your app, with a link to your Blood Glucose View Data screen. You'll also see Celebrations when you log blood glucose data for two or more days in a row, and when you come back and log after a break of seven days or more.

5.5 Adjust your target glucose range

PRECAUTION

Always consult with a physician before changing your target range.

The American Diabetes Association (ADA) recommends that you maintain a glucose range of 80-180 mg/dL.

You can also adjust your target glucose range.


1. Open the Onduo app.
2. Tap **Settings**  > **Target range** .
3. Tap **Your current range**, then scroll to change the maximum glucose reading for your target range.
4. Tap **Save**.

5.6 Review your blood glucose readings

To view a complete history of your day including events like exercise, steps, and medication, use the History tab:

1. Open the Onduo app.
2. Tap **History** .
3. Tap **Logbook**
4. Tap view Data within the Glucose category
5. Tap at the top of the screen to adjust your view by day, or week.
6. Tap the arrows at the top to see previous days or weeks.

On each day of your history, you can see a graph of your glucose readings, meals, and steps compared to that day's target range.

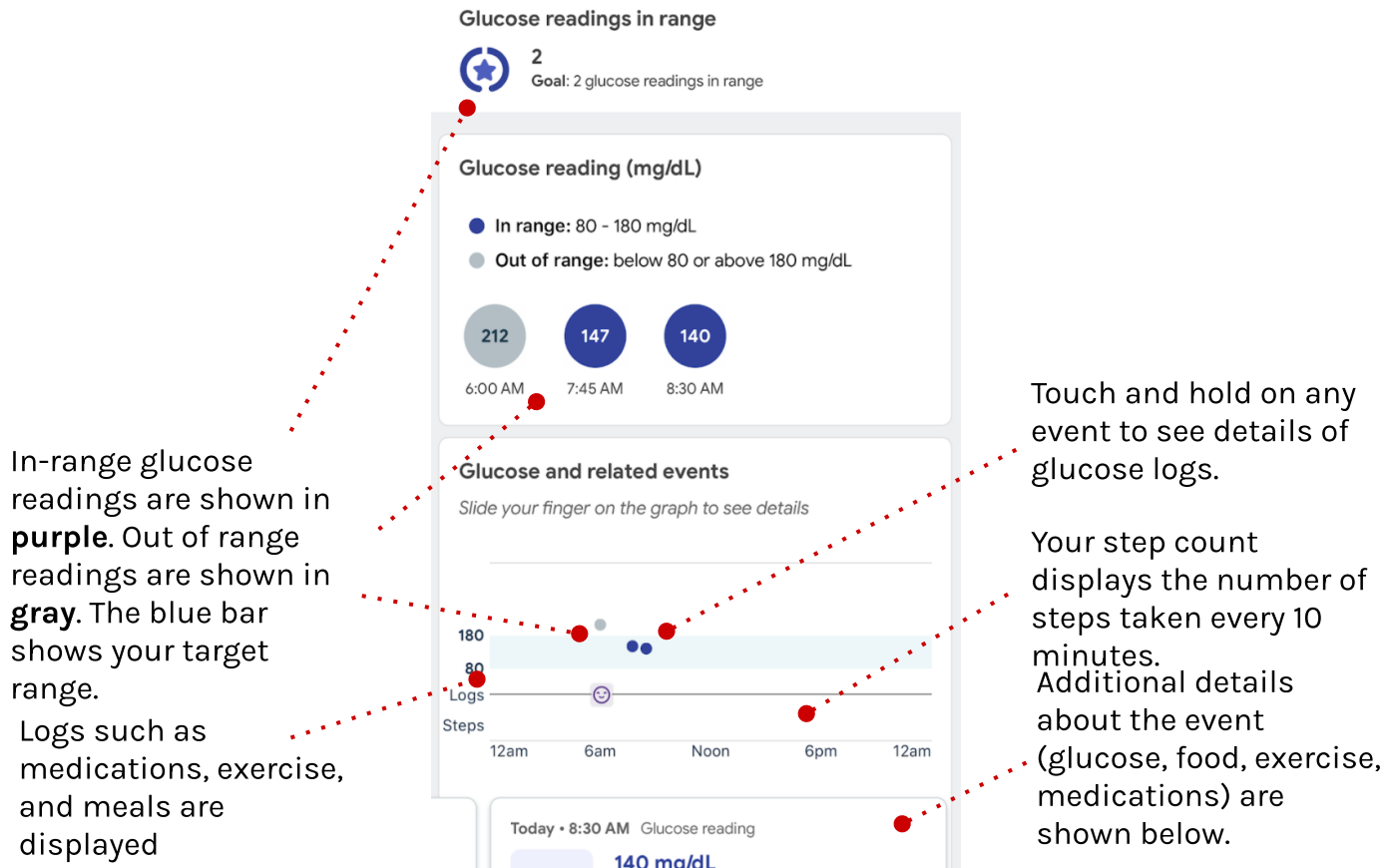You're not able to delete your history.

5.7 Manage your blood glucose meter (BGM)

If you receive, connect, and use a connected BGM as part of a program within Onduo, your glucose data becomes associated with your Onduo account. You are not able to disconnect a meter from your Onduo app yourself. To stop sending your blood glucose readings, call or email support at support@onduo.com or 1-833-HiOnduo (1-833-446-6386) and a support representative will help you disconnect the meter from your app. If you're participating in a clinical study, please contact the study support staff. Any readings that you logged up to the date you disconnect the meter will stay in your Onduo account.

Optionally, you can simply stop using your connected meter.

5.8 Blood glucose alerts

If you manually log a very low blood glucose reading, the Onduo app will alert you and provide advice on how to raise your blood glucose. For very low readings taken with your connected BGM, please pay attention to the alert displayed on the device. The Onduo app will not alert you for automatically saved readings.

Your Care Team may also check in with you via the Messages tab after a very low blood glucose reading is saved.

6.0 Continuous glucose monitor (CGM) (Type 2 diabetes only)

You may receive the opportunity to use a CGM as part of the Onduo Type 2 diabetes program.

PRECAUTION

Defer to your glucose device for primary glucose readings. If you believe the Onduo app is displaying incorrect glucose data, please refer to your glucose device.

For specific instructions for using the Dexcom CGM device, please refer to the CGM's instructions for use.

Note that glucose readings will be displayed on a 3-4 hour delay in the Onduo app.

6.1 What's a CGM?

A CGM utilizes a wearable device that inserts a sensor just under the skin to measure the amount of glucose in your interstitial fluid. Interstitial fluid refers to the fluid between the cells. Glucose levels in the interstitial fluid are similar to blood glucose levels, but the interstitial glucose levels may lag behind blood glucose levels by several minutes.

6.2 Connect your Dexcom CGM

With a CGM, a glucose reading is taken automatically every 5 minutes. When you connect Onduo to your CGM through the Apple Health app (iOS) or Dexcom API (Android), you can view and track your glucose readings at a three to four-hour delay without having to manually log this data.

To connect your CGM, follow all of the steps below.

6.2.1 Set up and insert your Dexcom G5 or G6 CGM

You will need the Onduo app and the Dexcom app to be able to receive your CGM data.

To set up and insert your CGM:

1. Download the Dexcom G5 or G6 app.
2. When the download is complete, open the Dexcom app.
3. Create an account and log in.
4. Use the instructions provided by Dexcom to set up and insert your CGM.

6.2.2 Connect your CGM data to the Onduo app

If you are part of a program in Onduo with a CGM, you can have glucose readings you take get saved and tracked in the Onduo app automatically using a cellular connection:

1. Open the Onduo app.
2. On the Home  tab, under the section titled "Connect your CGM," tap **CONNECT**.
If you don't see this, open the Onduo app **Settings** > **Connected meters** > DEXCOM_G5/G6 **CONNECT**.
3. Follow the on-screen instructions in the app to complete the connection process.

For specific instructions for using the Dexcom G5 or G6 device (such as insertion or on-device alerts), please refer to the Dexcom instructions for use.

6.3 Track and adjust your daily goal

While using a CGM, you should work towards a daily goal of keeping your glucose readings within your target range every day. You can see your progress on this goal in the **Day** section of the History tab.

6.3.1 Adjust your time in range daily goal and your target range

PRECAUTION

Always consult with a physician before changing your target range.

The American Diabetes Association (ADA) recommends that you maintain a glucose range of 80-180 mg/dL.

If it seems either too difficult or too easy to get any readings in your target range for the set amount of time, you can manually adjust your time in range.

1. Open the Onduo app.
2. Tap **Settings**  > **Time in range** .
3. Tap the time, then scroll to change your daily goal of number of hours in range.
4. Tap **Save**.

You can also adjust your target glucose range.

1. Open the Onduo app.
2. Tap **Settings**  > **Target glucose range** .
3. Tap **Current target range**, then scroll to change the maximum glucose reading for your target range.
4. Tap **Save**.

6.4 Review your glucose history

To view a complete history of your day including events like exercise, steps, and medication, use the History tab:

1. Open the Onduo app.
2. Tap **History** .

3. Tap **Logbook**
4. Tap view Data within the Glucose category
5. Tap at the top of the screen to adjust your view by day or week.
6. Tap the arrows at the top to see previous days or weeks.


See the progress

you've made on your
time in range goal for
the day.

In-range glucose readings
are shown in **purple**.

Out-of-range readings are
shown in **gray**. Readings
from your blood glucose
meter will be shown as
dots on your graph.

Logs such as
medications, exercise,
and meals are displayed

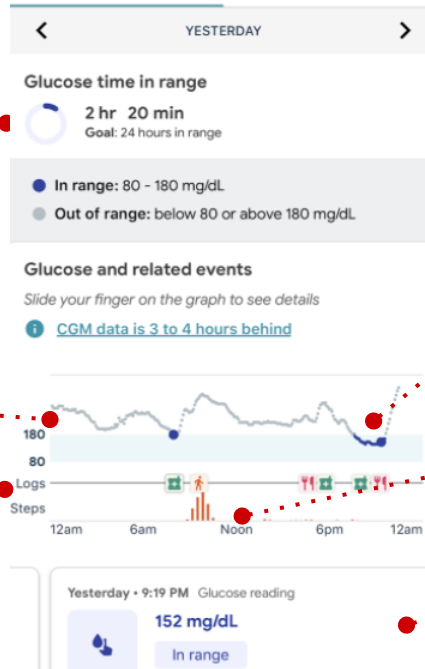Touch and hold on any
event to see details of
glucose data.

Your step count
displays the number of
steps taken every 10
minutes.

Additional details are
shown below the graph
or in these cards.

You're not able to delete your history.

6.5 Manage your Dexcom continuous glucose monitor (CGM)

If you receive and use a Dexcom CGM as part of a program within Onduo, your glucose data becomes associated with your Onduo app.

To disconnect your Dexcom CGM from the Onduo app:

1. Open the Onduo app.
2. Tap **Settings**  > **Connected meters** .
3. Tap **DISCONNECT** next to "Dexcom G6."

7.0 Home sleep test (HST)

Taking a home sleep test (HST) is a critical step in the process of getting an OSA diagnosis. If your physician prescribes you with a home sleep test after your evaluation video visit, the app will prompt you to enter your shipping information to order your HST. The test will be shipped to you. Follow the HST manufacturer's instructions for use.

Once you've completed your HST, your sleep physician will review your results. The Onduo app will prompt you to schedule a second video visit. You'll review your results and talk about your next steps with your sleep physician. If your sleep physician diagnoses you with OSA and prescribes PAP therapy, you're ready to move forward with the program.

8.0 PAP device

Your PAP device is an important part of your experience with the Onduo app if you're using the app to help manage obstructive sleep apnea. If you're enrolled in the sleep apnea program, you may receive a PAP Welcome Kit. Start by learning how to set up your PAP device, connecting it to the Onduo app, and reviewing information about your PAP use.

8.1 What's a PAP device?

PAP (positive airway pressure) therapy is a treatment for sleep apnea. When you use a PAP device, you wear a mask while you sleep. A tube connects the mask to a portable machine. The machine gently blows air into your airway. This air helps keep your airway open and supports normal breathing.

8.2 Set up and connect your PAP device

Wait until your PAP Therapy Kit arrives to begin setting up the PAP device and mask.

- PAP device set-up: The Onduo app will direct you to the manufacturer's electronic user guide to help you set up your PAP (for more details, refer to the manufacturer's IFU).
- Mask fit: Get tips for wearing the mask and easy breathing. Practice breathing with the mask, while sitting and lying down on a bed, and with your PAP device turned on.

8.3 Using the PAP device

The Onduo app will direct you to the manufacturer's electronic user guide to help you set up your PAP (for more details, refer to the manufacturer's IFU). Follow the instructions and activities in your app to help you adjust to PAP therapy, make PAP use a daily habit, and get comfortable falling asleep with your device. The Onduo app provides motivational tips to help you develop a habit of using PAP therapy and improve your sleep outcomes. Your device will sync your nightly PAP data with the app automatically - you don't need to upload anything manually after you use the device.

8.4 Review your PAP use

You can find detailed information on your PAP use in the **History section** of the Onduo app. This information is compiled into daily and weekly trends, and includes:

1. **Total PAP Usage:** This shows how much time, on average, you are using PAP during the night.
2. **Mask Leak:** Mask leak is an estimate of the total rate of air unintentionally escaping your mouth or mask. This can cause some discomfort, but does not mean the PAP isn't doing its job.
3. **AHI (apnea-hypopnea index):** The AHI reflects the average number of abnormal breathing events detected per hour of PAP use. Abnormal breathing events are scored when the device detects complete or partial stoppage of airflow through the nose and/or the mouth. Any events that you had while you were **not** using PAP are not included in this measure.

9.0 Medications

If your Care Team sets a medication logging goal, you can get helpful notifications on when to log your medications by adding them to the app.

PRECAUTION

This system is not intended to be used as a medication intake reminder system. Always defer to your prescribing physician's instructions. If you receive a medication logging notification at the wrong time, edit your medications list (refer to **Section 6.2**) and notify your Care Team.

Medication notifications remind you to log your medication intake. To turn on medication notifications, refer to **Section 9.1** to adjust your notification settings.

Follow your physician's advice if it conflicts with the app. The Onduo app is not intended as a substitute for your physician's guidance. Before making major lifestyle changes, please consult your physician.

If you take insulin, consult with your physician before making any lifestyle changes. You may need to modify the amount of insulin you take if you change your lifestyle, including food and exercise.

9.1 Add your medications

1. Open the Onduo app.
2. Tap **Settings**  > **Medications**  > **Add a new medication**.
3. Type in the name of the medication.
4. [Optional] Type in the dose of the medication.

5. Select which days you normally take the medication.
6. Tap **Add a time** to enter the times of day that you take the medication.
7. [Optional] Toggle logging reminders On or Off. By default, logging reminders are turned on.
8. Tap **Save**.

9.2 Edit and delete medications

1. Open the Onduo app.
2. Tap **Settings**  > **Medications** .
3. Tap the medication you want to edit or delete.
4. Adjust the name, dose, day, and time(s) you typically take the medication.
5. To remove times that you take the medication, tap .
6. Adjust logging reminders by toggling reminders On or Off.
7. Tap **Save**.
8. To delete the medication, tap **Delete Medication**.

9.3 Importing medications and medical supplies prescribed to you

The Onduo app may assist in adding medications and supplies for you by importing prescriptions that you fill. You will need to confirm that the imported prescriptions are correct.

1. Open the Onduo app.
2. On the Home  tab, under the section titled “New! Automate your med list,” tap **LEARN MORE**.

Note: This section is displayed only if you have not completed it in the past.

3. Follow the on-screen instructions in the app to opt in to the import process.

Note: The Onduo app may not be able to import all of your prescriptions. You have the option to manually add additional prescriptions (refer to **Section 9.1**).

9.4 Confirm medications imported from Surescripts or added by your Care Team

You will need to confirm that any medications imported, or added by your Care Team, are correct.

1. Open the Onduo app.

2. On the Home  tab, under the section titled “New prescription found!” tap **REVIEW NOW**.
3. Tap **ADD TO LIST** for the prescription you want to confirm.
4. Tap **YES** or **NO** if the new medication replaces an existing medication on your list.
5. Adjust the name, dose, day, and time you typically take the medication.
6. Tap **Save**.
7. If you do not want to add a new prescription to your list, tap **DON'T ADD**. Tap the reason for not adding the prescription.
8. You can also get to the new prescriptions by tapping **Settings**  > **Medications**  .

9.5 Log medications you have taken

Track when you took your medication, so you can see how it affects your glucose readings in your history.

1. Open the Onduo app.
2. Tap the task menu  > **Log medication**.
3. Tap the time to adjust when it was taken.
4. Check the box next to the medication or medications that you took.
5. Tap **Save**.

Note: You must add a medication (refer to **Section 9.1**) for you to be able to select and log it.

10.0 Food, exercise, weight, and mood logging

Log your meals, exercise, and mood to get a better understanding of how your health is impacted by your mood, activity and what you eat. This information is visible to both you and your Care Team—and you can work together to set goals for yourself.

10.1 Log food or a meal

1. Before you eat or drink, open the Onduo app.

2. Tap the task menu  **Log** > **Log meal**.

3. Tap the Photo banner  to add a photo.*

4. Type a description of what you ate and/or drank.*

5. Tap the time to adjust when you had this meal.

6. [Optional] Tap either **Breakfast**, **Lunch** or **Dinner** to specify if it was one of the three major meals.

7. [Optional] Tap **Choose a location** and pick either **Find on map** or **At home** depending on where you had the meal.

8. Tap **Save**.

9. The first time you log your food, you'll see a "Celebration" message in your app, with a link to your Food View Data screen. You'll also see Celebrations when you log your food for two or more days in a row, and when you come back and log after a break of seven days or more.

**Either a photo or a food description is required to save this log. Both are not required.*

Once logged, it's not possible to edit or delete logged meals.

10.2 Log exercise

You can manually log different types of exercise:

- 
-
- 
- >
- Log exercise.**
1. Tap the task menu  > **Log exercise.**
 2. Tap the icon to choose the type of exercise you completed.
 3. Tap the numbers to adjust your start time and duration of exercise.
 4. Tap **Low** to change the intensity of the exercise.
 5. [Optional] Include a note about the activity for yourself and your Care Team.
 6. Tap **Save**.
 7. The first time you log your exercise, you'll see a "Celebration" message in your app, with a link to your Exercise View Data screen. You'll also see Celebrations when you log your exercise for two or more days in a row, and when you come back and log after a break of seven days or more.

You aren't able to connect your own activity tracker (like a FitBit). You can track steps (refer to **Section 11.0**).

10.3 Log mood

You can manually log how you are feeling:

- 
-
- 
- >
- Log mood.**
1. Tap the task menu  > **Log mood.**
 2. Tap the icon that best describes how you are feeling.
 3. Write about a few things that happened this week, and describe how they made you feel.
 4. Tap the numbers to adjust your start time and date.
 5. Tap **Save**.
 6. The first time you log your mood, you'll see a "Celebration" message in your app, with a link to your Mood View Data screen. You'll also see Celebrations when you log your mood for two or more days in a row, and when you come back and log after a break of seven days or more.

10.4 Log weight

You can manually log your weight:

- 
1. Tap the task menu **Log** > **Log weight**.
 2. Enter your weight in lbs.
 3. Tap the numbers to adjust date and time.
 4. Tap **Save**.

10.5 View your logged events

You can view logged events such as meals you've logged, steps you've taken, and exercise you did.

1. Open the Onduo app.
2. Tap **History** .
3. Tap **Logbook** and select **View Data** in the **Food/drink, Movement, Medication, or Mood** category to see your logged events in more detail. (If you are in the Onduo Mental and Behavioral Health program, the option will be titled "Mental Health" instead of Mood.)
4. Scroll left and right through the cards on the bottom of the page to see information about each day.
5. Use the arrows at the top to see logged events from the past.

10.6 Food logging suggestion

When you are wearing a CGM the software application may detect that a meal has been consumed (by examining your CGM data) and suggest to log the meal. In order to log the meal based:

1. Tap on the notification to open the Onduo app.
2. Tap **Yes, let's log it**
3. Enter a photo, description, adjust time/date and other relevant details
4. Tap **Save**

11.0 Step tracking

Track your steps in the Onduo app to better understand how your movement impacts your glucose readings. This information is visible to both you and your Care Team—and you can work together to set goals for yourself. **Please note:** step tracking has been temporarily disabled for Android users.

11.1 Track steps and adjust your step goal

Your steps are automatically tracked by data gathered from your phone. You can change the number of steps to aim for per day:

Your phone automatically tracks the number of steps you take each day. To share this data with Onduo, allow the Onduo app to get step data from Google Fit or Apple Health:

1. Open the Onduo app.
2. Tap **Settings**  > **Steps** .
3. Follow the on-screen instructions to get steps data from Google Fit or Apple Health.
If you use Android, you may need to download Google Fit.

You can change the number of steps to aim for per day:

1. Open the Onduo app.
2. Tap **Settings**  > **Steps** .
3. Move the slider to change your daily step target or tap the number to manually adjust.
4. Tap **Save**.

11.2 View your steps

You can view steps you've taken.

1. Open the Onduo app.
2. Tap **History**  > **Logbook**.
3. Tap on **View Data** in the **Movement** category.

4. Tap at the top of the screen to adjust your view by day or week.
5. Tap the arrows at the top to see previous days or weeks.

Steps are also shown on the Home tab so you can track your progress for the day.

12.0 Your Journey: Activity sets, activities, checkpoints, insights, and education

From your Home tab, you can choose a journey that's important to you. Journeys are designed to help you achieve your health goals by having you work towards your ideal glucose range, hitting your target blood pressure level, or improving your mental and behavioral health. Each journey is made up of activity sets. Each activity set contains up to eight activities for you to complete.

12.1 Activity Sets

Activity sets consist of activities, informational content, and challenges designed to help you develop healthy habits for managing Type 1 diabetes, Type 2 diabetes, hypertension, or your mental and behavioral health. Activity sets will appear in the "Your program" section on the Home tab. Each activity set lasts about a week. You may be able to see what's coming next or view your progress anytime by clicking "View full journey."

12.2 Learn about activities

Activities focus on different aspects of your health to help you improve by tracking actionable tasks. For example, you could have an activity around logging your glucose reading before and after a meal to see how your food affects your glucose readings. Alternatively you might learn about small ways to change your evening routine to help you get a better night's sleep.

Each activity lasts about a week, but your daily goal will always be there for you. Some activities can make achieving your daily goal easier!

12.3 Learn about insights and educational materials

Your Care Team may surface insights and educational materials to you in the app. The more information the Care Team has about your habits and health, the better they can provide insights and education to you.

These insights might include basic information about diabetes, high blood pressure, or mental and behavioral health. Or you might get personalized information about how your glucose or blood pressure readings are affected by your breakfast.

12.4 Reminders for activities

If you have notifications turned on, you'll be reminded what tasks are needed to complete the activity. The app will provide you with updates on the progress of your activity on the Home tab.

12.5 Checkpoints

Checkpoints help track your progress on a clinical measure connected to a clinical condition. Checkpoints typically appear in the “For you” section. Each checkpoint has a set duration which is generally around a week.

13.0 Video visits

You may be given the opportunity to meet with Onduo Care Team members through a video visit.

13.1 Schedule a visit

Your Care Team will let you know when it's time to schedule a video visit.

1. Open the Onduo app.
2. Navigate to the message from your coach which shares a link for you to schedule an appointment.
3. Click on the link and follow the instructions displayed on-screen to schedule, modify, and join your appointment.

13.2 (Old flow - No longer in use) Schedule a visit

This flow will no longer be used going forward, but instructions are provided for reference for members who had it triggered previously and have not yet completed it.

Your Care Team will let you know when it's time to schedule a video visit.

4. Open the Onduo app.
5. Navigate to the message from your coach which shares a link for you to schedule an appointment.
6. Click on the link to schedule your appointment.
7. Click the "Find a time" button on the scheduling card on the home screen. Tap "Show more" to view available appointment times over the next three months.
8. Choose an appointment time that works for you.
9. Indicate if you would like a Spanish interpreter to join your video visit.

13.3 (Old flow - No longer in use) View and modify your appointments

This flow will no longer be used going forward, but instructions are provided for reference for members who had it triggered previously and have not yet completed it.

See section 18.4 My Appointments

13.4 (Old flow - No longer in use) Join a video visit

This flow will no longer be used going forward, but instructions are provided for reference for members who had it triggered previously and have not yet completed it.

1. Open the Onduo app.
2. Click the "Check in" or "Join now" button at the top of the home page. This button will be available from 10 minutes before the appointment starts until 5 minutes after the appointment ends, or until you exit the visit.
3. Follow the prompts if shown to set up your camera and microphone.
4. Click the "Enter visit" button to join the video visit.

14.0 Blood pressure

If you are enrolled in the hypertension management program, you may be given the opportunity to track your blood pressure in the Onduo app. We recommend you use a blood pressure monitor that connects to the Onduo app. Please view the list of devices in Settings > Connected meter. You can also choose to log your blood pressure and pulse readings manually.

Start by learning how to connect your monitor and review your readings.

PRECAUTION

Defer to your blood pressure monitor for primary blood pressure and pulse readings. If you believe the Onduo app is displaying incorrect blood pressure or pulse data, please defer to the readings on your blood pressure monitor. Please contact the Onduo support team if you believe the app is displaying incorrect data. For specific instructions for using your blood pressure monitor, please refer to the monitor's instructions for use.

14.1 Indications for use/Target population

The Onduo Hypertension Management Feature is intended to be used by individuals who:

- Are aged 18 and older,
- Have been diagnosed with Hypertension, and
- Are using iOS and Android mobile operating systems.

The Onduo Care Team Console is intended to be used by healthcare professionals and care team users.

14.2 What's a blood pressure monitor?

You may receive a blood pressure monitor as part of Onduo's hypertension management program. Or you may be able to connect a blood pressure monitor that you already own. These devices measure how hard your blood pushes against the walls of your blood vessels as it flows.

14.3 Connect your blood pressure monitor

You can connect most blood pressure monitors to the Onduo app to automatically save and track your blood pressure and pulse readings in the Onduo app.

1. Open the Onduo app.
2. On the Home  tab, under "Pair your blood pressure monitor," tap **Pair**.
If you don't see this, open the Onduo app **Settings** > **My Devices** > Add a device.
3. Choose your device from the list of approved BP devices.
4. Follow the instructions in the app, and remember to insert batteries into your monitor.
5. The app will confirm once your monitor has been successfully connected. If the connection was unsuccessful, you'll see a failure message with a Retry button.
6. The first time a blood pressure reading is received from your monitor, you'll receive a "Celebration" push notification with a link to your Blood Pressure View Data screen. You'll also see Celebrations when your app receives blood pressure data for two or more days in a row, and when data is received after a break of seven days or more.

If you have questions, or if you don't see the option to "Pair your blood pressure monitor," please contact your Care Team. Or, try opening the Onduo app **Settings** > **My Devices** > Add a device.

14.4 Manually log a blood pressure reading

You can log blood pressure readings manually:

1. Open the Onduo app.
2. Tap the task menu  > **Blood pressure**.
3. Type in your systolic and diastolic readings (recorded as mmHg).


4. You can optionally enter your pulse rate (recorded as BPM, beats per minute), if you know it.
5. Adjust the date and time of the reading, if needed.
6. Tap **Save**.
7. The first time you manually log a blood pressure reading, you'll see a "Celebration" message in your app, with a link to your Blood Pressure View Data screen. You'll also see Celebrations when you log blood pressure data for two or more days in a row, and when you come back and log your blood pressure after a break of seven days or more.

14.5 Review your blood pressure and pulse readings

To view a complete history of your day including events like exercise, steps, and medication, use the History tab:

1. Open the Onduo app.
2. Tap **History**  **> Logbook**.
3. Tap **View Data** in the **Blood Pressure** category.
4. Tap at the top of the screen to adjust your view by day or week or month.
5. Tap the arrows at the top to see previous days or weeks or months.

On each day or week of your history, you can see a table of your blood pressure and pulse readings. In the weekly and monthly views, you can also see your average blood pressure for that week or month.

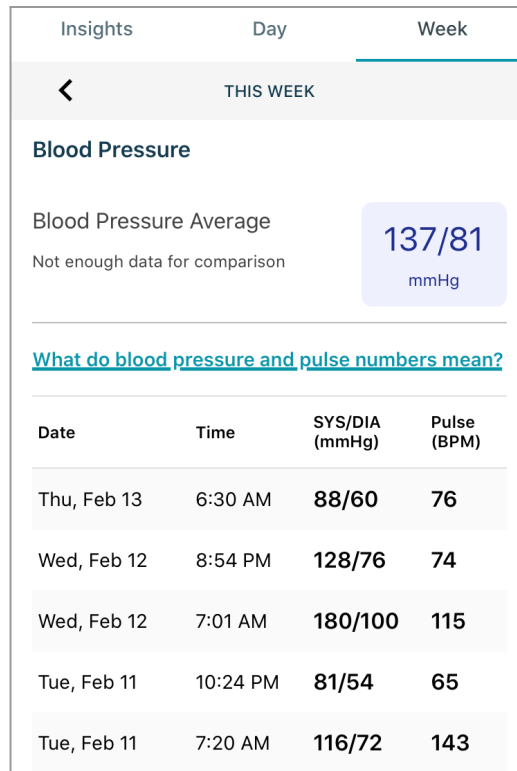You're not able to delete your history.

14.6 Finding your baseline blood pressure

You may be given the opportunity to participate in a blood pressure checkpoint by your Care Team. Checkpoints show you how your blood pressure is doing and help you see if you might want to make changes.

If you opt-in to the blood pressure checkpoint, you'll be asked to check your blood pressure in the morning and evening for 3 days. You can keep track of your checkpoint progress through a card on your home screen.

Once you complete a checkpoint, you will be able to see your average blood pressure in the morning and evening, and your overall average blood pressure. Your care lead will review this with you and talk about a blood pressure plan.

14.7 Manage your blood pressure monitor

If you use a blood pressure monitor as part of a program within Onduo, your blood pressure and pulse data becomes associated with your Onduo app.

To disconnect a blood pressure monitor from the Onduo app:

1. Open the Onduo app.
2. Tap **Settings**  > **My Devices**  .
3. Tap the device you wish to disconnect, and tap **Forget device**.

14.8 Blood pressure and pulse reading alerts

Very low or high blood pressure or pulse readings, whether read directly from a blood pressure monitor or manually logged, will trigger alert notifications. The alert will prompt you to take another reading after resting for 5 minutes. If the reading is still very low or high, the alert will prompt you to talk to your care team. Your Care Team may also check in with you via the Messages tab when this occurs.

15.0 Obstructive Sleep Apnea

The Onduo Obstructive Sleep Apnea (OSA) program is designed to help people adjust to positive airway pressure (PAP) therapy, and to build and maintain healthy sleep behaviors and routines.

The OSA program can help you:

- **Get an OSA diagnosis:** Complete a sleep questionnaire, video visits with a sleep physician, and a home sleep test to help determine whether you have an official OSA diagnosis.
- **Set goals:** Think about how OSA affects your life, and what you'd like to change the most.
- **Get your mask recommendation:** Answer a few questions and get your recommended PAP mask style, based on your needs and preferences.
- **Get your PAP welcome kit:** Order your full PAP welcome kit by following our step-by-step guide.
- **Set up your PAP therapy:** The Onduo app will direct you to the manufacturer's electronic user guide to help you set up your PAP (for more details, refer to the manufacturer's IFU). Trained Onduo coaches and support staff can help with any PAP equipment or adjustment issues.
- **Feel supported:** Work with a Care Lead to make PAP therapy part of your daily routine.
- **Track progress:** View insights about your PAP usage to help you stay on track.

15.1 Indications for Use/Target population

The Onduo Sleep Apnea Program is intended to provide an Obstructive Sleep Apnea (OSA) management program. Members enrolled in the OSA program will be assessed for obstructive sleep apnea through home testing and professional evaluations. After receiving an official OSA diagnosis, if PAP therapy is recommended, members gain access to PAP equipment and instructions for set up and adjustment, care team support, and guidance and education.

The OSA program will offer support and strategies to help manage diagnosed members' symptoms through self guided digital activities and lifestyle management skill building. These activities are provided via in-app education on OSA therapy, as recommended by their Care

Leads. Members can message Health Coaches and get connected with qualified healthcare professionals to engage in clinical consultation via telehealth visits.

The Onduo Sleep Apnea program is intended to be used by individuals who are:

- Aged 18 and older,
- Using iOS and Android mobile operating systems,
- Experiencing symptoms of OSA, and
- Not currently under treatment for OSA.

15.2 Is Onduo's Sleep Apnea program right for me?

Onduo's Sleep Apnea program is not right for you if you:

- Use supplemental oxygen,
- Have a history of stroke,
- Use opioid medications,
- Have severe insomnia,
- Have central apnea,
- Have congestive heart failure,
- Have neuromuscular disease, or
- Have moderate to severe pulmonary disease.


Onduo Sleep Apnea program members must participate in a home sleep test and use PAP therapy as part of the program. If your sleep physician decides that a home sleep test or PAP therapy is not the best fit for you, you will not be enrolled in the Onduo Sleep Apnea program.

15.3 Getting an OSA diagnosis

Members of Onduo's Sleep Apnea program must receive an OSA diagnosis through the program and be prescribed PAP therapy by the sleep physician as the appropriate next step before beginning PAP therapy. Follow these steps to get a diagnosis and get started:

1. **Confirm eligibility:** Your initial Onduo outreach materials contain a link to the Onduo website. Follow the steps to confirm your eligibility and coverage. Next, you'll be prompted to create a log-in and download the Onduo app.

2. **Complete onboarding in the app:** Once you've downloaded and logged into the app, you'll answer some onboarding questions and get matched with a Care Lead.
3. **Pre-diagnosis:** Next, you will receive an introduction to the OSA program through the Onduo app. You'll be asked some screening questions in order to continue through the diagnosis process.
4. **Diagnosis:** There are four steps to take in order to get an OSA diagnosis:

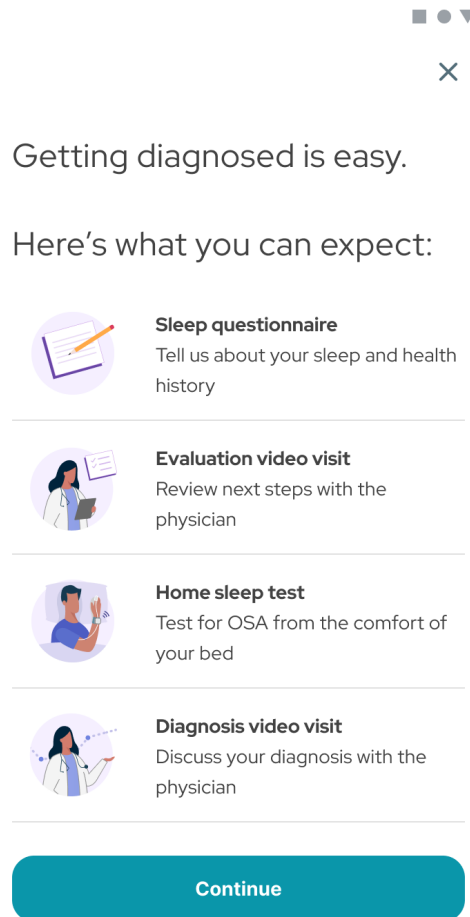Diagnosis screen for OSA program

- a. **Sleep questionnaire:** You'll answer questions about your sleep and health
- b. **Telemed visit #1: Evaluation:** You'll schedule a video visit with a sleep physician. The visit will happen through the app (don't worry, you'll get a reminder!).

Select a date and time ✕

Dr. Ava Nilsson

Program physician

Select a different doctor ▾

Choose from the times below (PT):

Wed, Aug 26, 2020

9:00 AM	9:30 AM	10:00 AM
11:30 AM	12:00 PM	12:30 PM
1:00 PM	2:00 PM	2:30 PM

Thu, Aug 27, 2020

11:30 AM	12:00 PM	12:30 PM
1:00 PM	2:00 PM	2:30 PM

Telemed scheduling screen for OSA program

Video visit screens for OSA program

- c. **Home sleep test:** If your sleep physician decides a home sleep test (HST) is right for you, the app will prompt you to order a home sleep test. The test will be shipped to you. Follow the test manufacturer's directions for setting up and using this device.
- d. **Telemed visit #2: Diagnosis:** Once your sleep physician has reviewed your HST results, the app will prompt you to schedule a second video visit. You'll review your results and talk about your next steps. If your sleep physician diagnoses you with OSA and prescribes PAP therapy, you're ready to move forward with the program.


15.4 Welcome Kit and PAP Therapy set up

Once you have received your formal OSA diagnosis and PAP therapy prescription from your sleep physician, you can get started with PAP therapy. Follow these steps to get set up:

1. **Build your PAP Welcome Kit:** Your app will prompt you to build your PAP Welcome Kit.


You'll answer questions about your goals for PAP therapy and choose the right equipment

for your kit. The kit will include materials for PAP therapy, including mask cushions, air filters, mask frame, tubing, headgear, PAP device, and other materials. Your Welcome Kit will be shipped to you via postal mail.

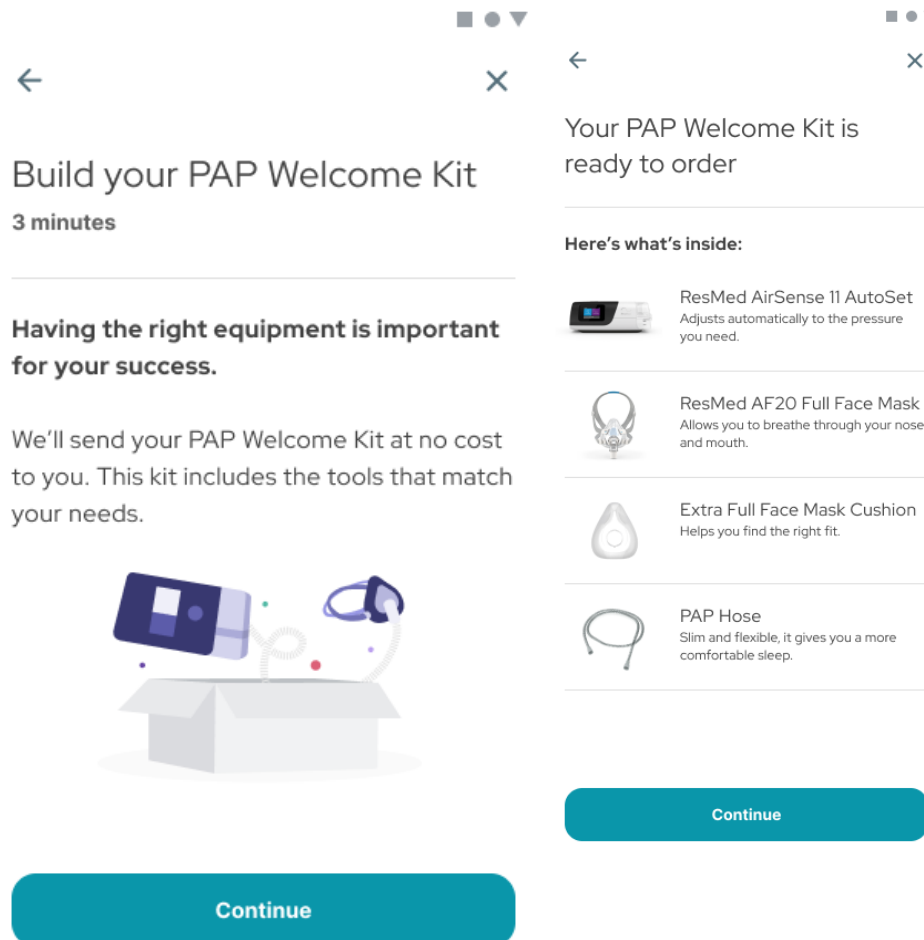PAP Welcome Kit screens for OSA program

2. **Learn while you wait:** While you wait for your Welcome Kit to arrive, you can learn more about OSA and PAP therapy in your app.
3. **Initial PAP setup:** When your Welcome Kit arrives, the Onduo app will direct you to the manufacturer's electronic user guide to help you set up your PAP (for more details, refer to the manufacturer's IFU).. You'll learn how to operate your device and fit your mask, and you'll get guidance on your first night using PAP therapy.

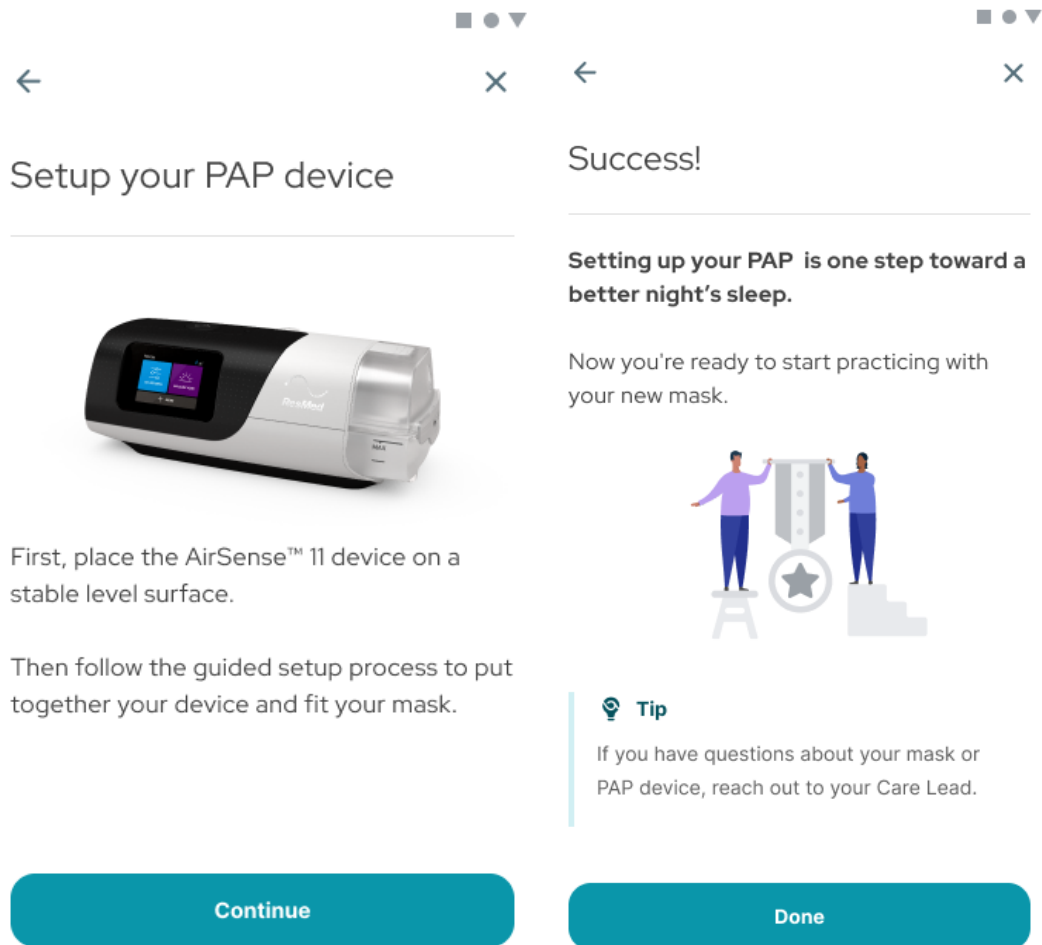Initial PAP set up screens for OSA program

4. **First week with PAP therapy:** Once you've set up your PAP device, practiced with device and mask, and experienced a first night using PAP, you'll get seven consecutive nights of guided audio activities to help you relax and use PAP effectively. You'll be prompted to reflect on each activity in the morning. At the end of the week, you'll see a summary of your responses along with your PAP usage data.

15.5 Using PAP therapy

PAP therapy is a critical part of the OSA program. PAP therapy involves wearing a face or nasal mask during sleep. The mask is connected to a hose, which is connected to a pump on the PAP device. The device pumps filtered air through the mask to keep your airway open while you sleep.

It's important to learn how to use your PAP device comfortably and confidently, and use it every night. Follow the instructions and activities in your OSA program to guide your breathing and get comfortable falling asleep with PAP. You'll also find motivational tips to help you develop a habit of relaxing with and using PAP therapy and improve your sleep.


15.5.1 Schedule first night and first week activities

Once you have scheduled your first night and/or first week activity, you'll automatically receive mobile notification reminders through your app. These notifications will prompt you to complete activities that support your PAP use.

15.5.2 PAP therapy tracking

Your PAP device works with the Onduo OSA Program to track information about your PAP use. You can find daily and weekly trends under PAP Therapy in the **History** section of the app. This information includes:

- **Total PAP Usage:** PAP Usage shows much time, on average, you're using PAP during the night.
- **Mask Leak:** Mask leak is an estimate of the total rate of air unintentionally escaping your mouth or mask. This can cause some discomfort, but it does not mean the PAP is not doing its job.
- **AHI:** The AHI (apnea-hypopnea index) on the PAP device reflects the average number of abnormal breathing events detected per hour of use. Abnormal breathing events are scored when the PAP device detects partial or complete stoppage of airflow through the nose and/or mouth. Any events you might have while not using PAP are not included in this measure.

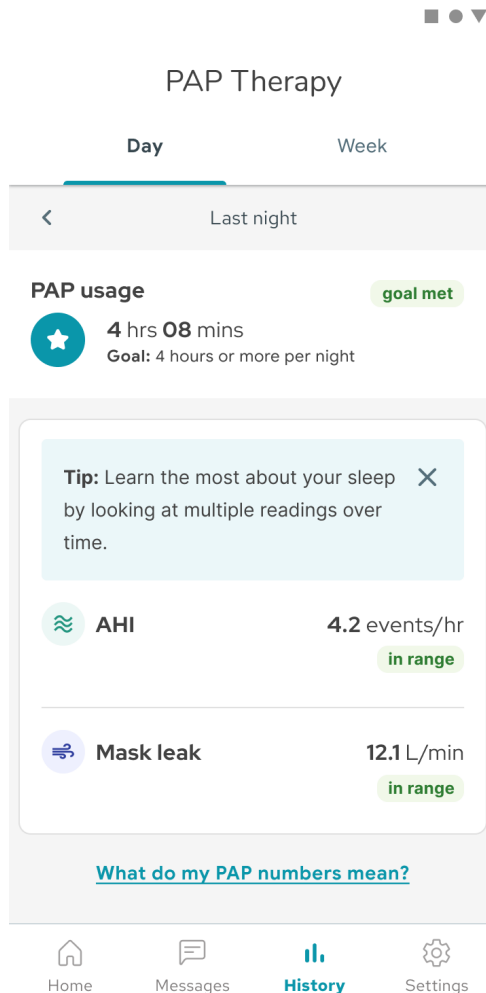PAP therapy history screen for OSA program

15.6 Using the Onduo app

The Onduo app has four key navigation tabs:

- **Home:** Access all OSA program activities, including learning about OSA and PAP therapy, PAP device and mask setup, guided relaxation, and healthy sleep routines.
- **Messages:** Send and receive messages from a Care Lead for personalized support through the journey with PAP therapy.
- **History:** Track PAP therapy data and insights, including PAP usage, apnea-hypopnea index (AHI), and Mask Leak.
- **Settings:** View your care team and your connected PAP device.

15.7 Your guided learning activities

Each activity in the OSA program includes a learning module, review questions, and reflection.

These activities can help you to:

- Understand OSA, and its effects on your health and quality of life,
- Learn about PAP therapy, and
- Get tips on day-to-day activities, such as your bedtime routine.

15.8 Messaging with your Care Team

If you are having any difficulties with your PAP device or PAP therapy, you can message your coach for assistance. You can also view your care team information under Settings > My care team to reach your provider and/or Durable Medical Equipment (DME) group.

15.9 Ordering replacement parts for your PAP device

The components of your PAP equipment can wear out, which can make your therapy less effective. You can use the app to learn about resupply and to contact your DME.

16.0 Mental and Behavioral Health

The Onduo Mental Behavioral Health program is intended to screen members for mental health conditions and facilitate appropriate interventions to manage their symptoms. Enrolled members will be provided with tools to manage their anxiety, depression or related symptoms through self guided digital activities and lifestyle management skills via in- app education as recommended by their care leads. As needed, members will be connected with recommended qualified healthcare professionals most appropriate to their needs to engage in mental health clinical consultation via telehealth visits.

16.1 Indications for use/Target population

The Onduo Mental Health Management feature is intended to be used by targeted toward individuals who are:

- Aged 18 and older,
- Using iOS and Android mobile operating systems, and
- Experiencing a mental health need, with or without an existing diagnosis, or
- Onduo users currently enrolled in order to improve outcomes in another condition, such as Diabetes or Hypertension.

The Onduo Mental Health Management feature is not intended for individuals suffering from high risk mental health disorders requiring emergency care, including recent suicide attempt and/or psychiatric admission, within 3 months.

16.2 Is Onduo's Mental and Behavioral Health right for me?


Onduo's Mental and Behavioral Health (MBH) program is meant for people who want to improve their mental health and well-being. You don't need to have an existing diagnosis. But Onduo's MBH program won't be right for everyone, and is **not** right for you if:

- You're suffering from symptoms of depression or another mental health disorder that can't be managed in Onduo and require a higher level of care (as determined by an Onduo licensed mental health care provider);
- You require emergency care for your mental health;
- You've been admitted to a hospital for psychiatric reasons in the past 3 months; or

- You’ve attempted suicide in the past 3 months;
- You live in a U.S. state where Onduo’s Mental and Behavioral Health program is not currently supported;

16.3 Getting started


Your coach will send you a card to opt in to the MBH program. This will also show a push notification on your phone. Tap the “Let’s talk about mental health” card on the homescreen of the app to get started.

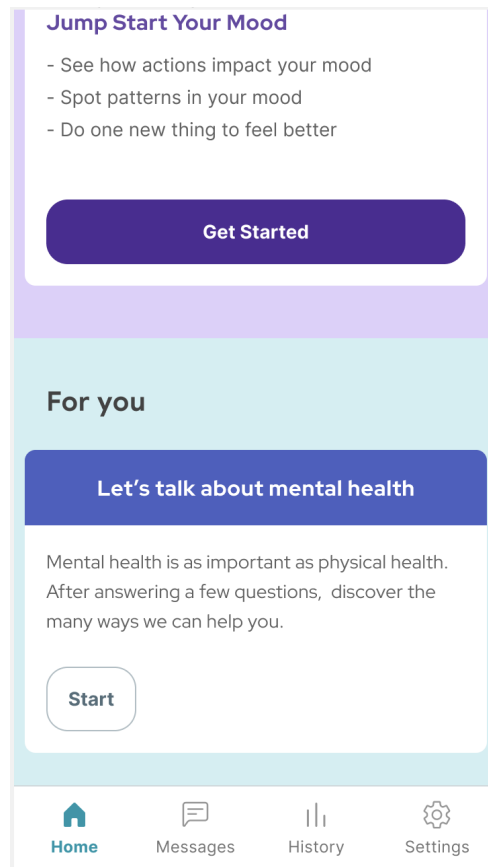Home Screen's Opt-in card for MBH program

First, we'll give you an overview of the program, then ask you to confirm you'd like mental health support from Onduo. Next, we'll ask a quick question about your medical history as it pertains to mental health, then ask you to provide an emergency contact, to make sure you're supported and safe. (You can always edit this emergency contact later in Settings).

After that, you'll answer the **PHQ-9 survey**. This survey helps us to evaluate your depression symptom severity. It will ask you nine questions about how you've been feeling lately.

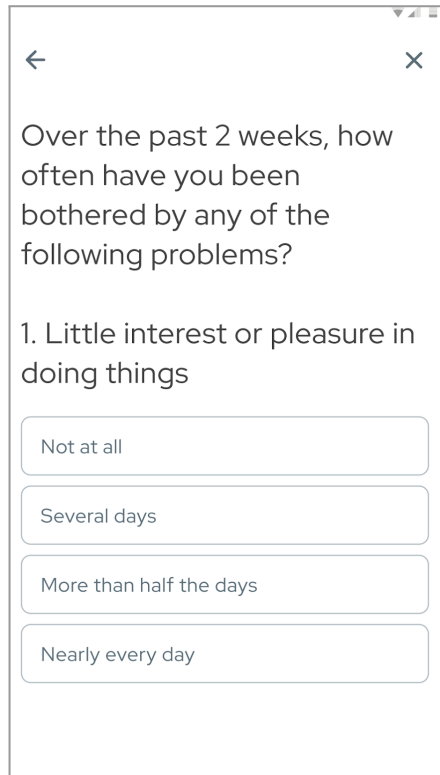← ×

Over the past 2 weeks, how often have you been bothered by any of the following problems?

1. Little interest or pleasure in doing things

Not at all


Several days

More than half the days

Nearly every day


PHQ-9 survey question screen

Your answers to the PHQ-9 will be used to create your **depression score**. Your depression score is a number between 0 and 27. Your score tells a licensed mental health care provider if you have symptoms of depression. Symptoms can be mild, moderate, or severe.

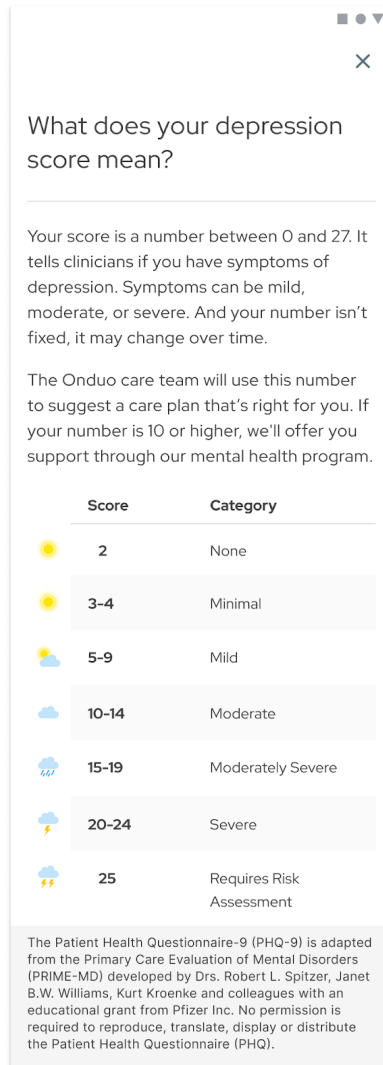Depression score screen

Your depression score can change over time. In the History section of the app under “Mental health,” you can compare your current score against past scores to track how you’re doing.

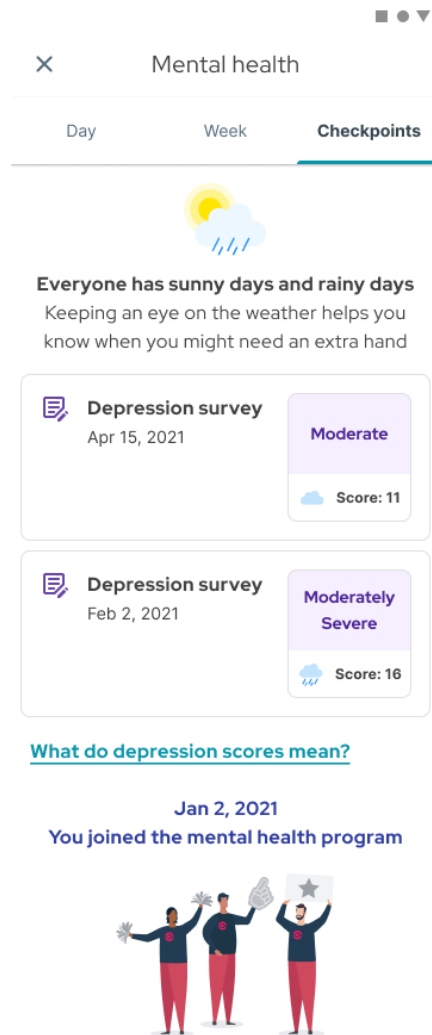
Depression checkpoint screen

The PHQ-9 is not intended to provide a diagnosis of depression, or to replace the judgment of a healthcare professional. The PHQ-9 is just one way to screen for depression but does not screen for other mental health conditions. Your depression score may become part of your health record. It may be used by a coach or a licensed mental health care provider to support and inform their evaluation and treatment of you. You should talk to your doctor if you have any questions about the PHQ-9 or your depression score.

Once you've finished the initial PHQ-9 survey questions, you'll be able to start your journey.

16.4 Your journey


As discussed in Section 9.0, journeys are designed to help you build new skills and habits, and achieve your health goals. Each journey is made up of activity sets. Each activity set contains up to eight activities for you to work on, and takes about a week to complete. You can choose which activities are right for you.

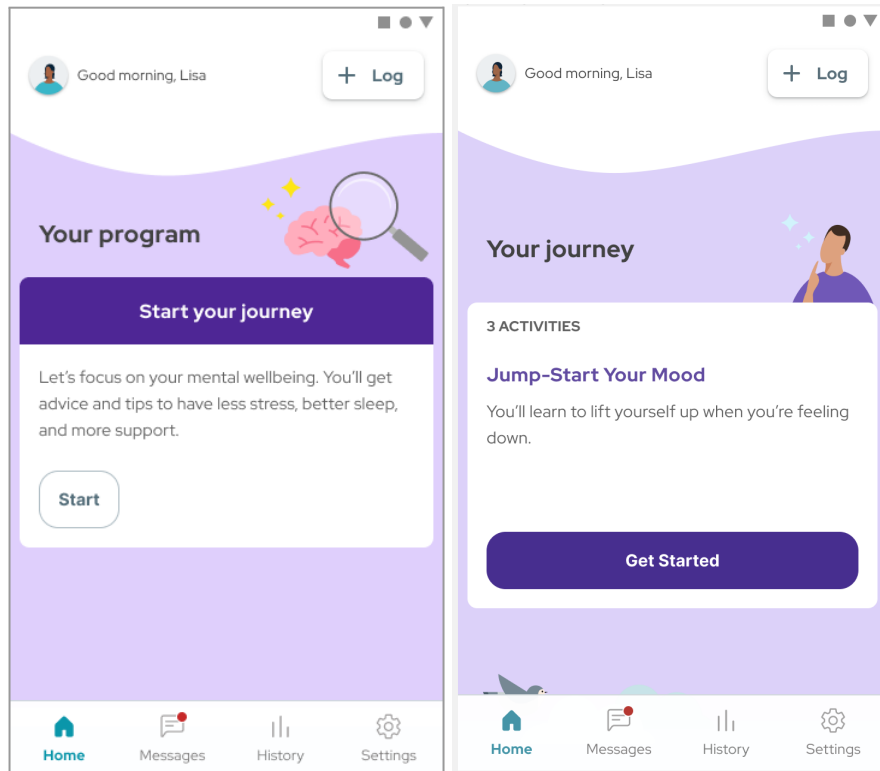Start your journey screen

16.5 Messaging with your Care Team


You have access to an Onduo Care Team through your program. You can send your Care Team a message anytime by clicking on the Messages button at the bottom of the screen. If you have an unread message, you'll see a red dot on the Messages icon. Your Care Team may send you content, articles, or suggestions. They may also check in to see how you're doing and feeling.

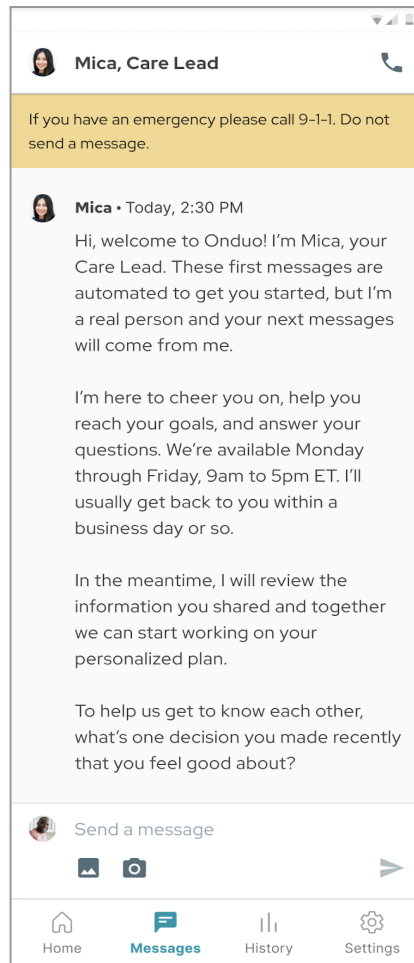Messages screen

16.6 Mental Health Checkpoints

From time to time, we'll check in and ask you to answer the PHQ-9 survey again. This helps inform your Care Team about how you're doing, and whether your program needs any changes to keep you moving forward. (For details on the PHQ-9 survey, see the "Getting Started" section above.)

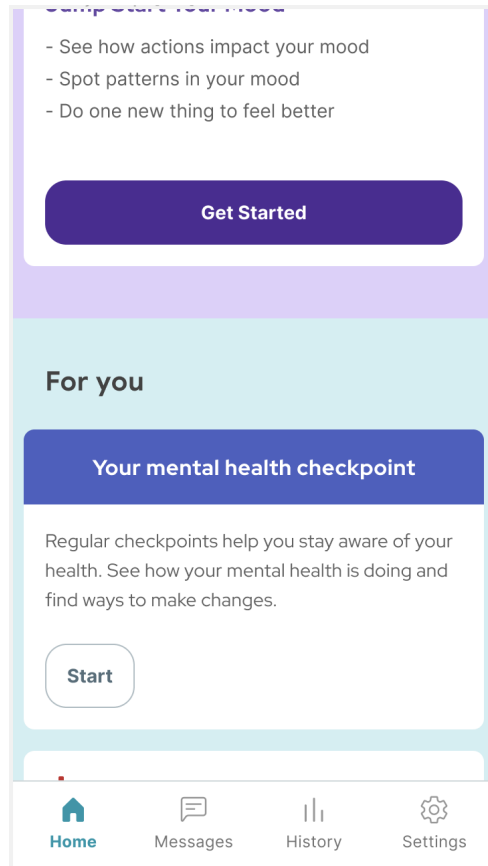Mental health checkpoint screen

16.7 Video visits with clinicians

You may get the opportunity to have a video consultation with a licensed mental health care provider as part of your program. During this consultation, you'll talk about setting up further video appointments with a licensed mental health care provider. Your consultation can help you decide which type of provider is best for you. Once you choose, you'll receive information on your home screen about how to set up your first appointment.

17.0 Type 1 Diabetes

The Onduo Type 1 Diabetes program is intended to support members with Type 1 diabetes with help from their healthcare professionals and wellness/care teams. The program provides review, analysis and evaluation of glucose data from compatible FDA-cleared blood glucose meters and FDA-approved continuous glucose monitoring systems in support of an effective diabetes management program.

17.1 Intended use

The Onduo Diabetes Management Software Platform is intended to provide individuals with Type 1 Diabetes, their healthcare professionals and wellness/care teams, a review, analysis and evaluation of retrospective glucose data from compatible FDA-cleared blood glucose meters and FDA-approved continuous glucose monitoring systems in support of an effective diabetes management program.

The Onduo Diabetes Management Software Platform is not intended to provide treatment decisions nor is it to be used as a substitute for professional healthcare advice. It is designed to help Onduo members manage their diabetes, symptoms and treatment of their health conditions through self-guided activities connecting with a CDCES or Diabetes Coach who are certified healthcare practitioners trained in leading care for insulin dependent diabetes members.

17.2 Target population

The Onduo Type 1 Diabetes program is intended to be used by individuals who:

- Have insulin-dependent diabetes (Type 1), and
- Are aged 18 and older, and
- Are using iOS and Android mobile operating systems.

The Type 1 Diabetes program is not intended for use by

The Onduo Type 1 diabetes management program is not intended for use by individuals who:

- have cystic fibrosis

- have cirrhosis or liver failure
- have Alzheimers
- are currently pregnant
- are institutionalized
- are receiving hospice care
- are receiving active chemotherapy treatment for cancer

17.3 Getting started

You'll first register for Onduo and download the Onduo app. You will take a survey about your health, needs, and goals. This information will help us match you with a Care Specialist. Once you're set up in the app, you can set up a phone and/or telemedicine visit with your Care Specialist to talk about your diabetes management. You will receive ongoing coaching from your care team. You'll also receive personalized diabetes education content that's tailored to your needs and goals.

17.4 Devices and supplies

Depending on your needs and preferences, you may receive the following supplies:

- **Blood glucose meter (BGM):** BioTel Blood Glucose Meter
- **Continuous glucose meter (CGM):** Your Onduo physician may prescribe a CGM device, as needed.
- **A1C kits:** You may receive A1C kits every 3-6 months, or as needed, to collect your A1C results.

17.5 Video visits with Onduo PC

A clinical pharmacist may refer you to the Onduo PC Physician for a telemedicine consultation. You must consent to being seen by the Onduo PC before an appointment can be scheduled. You also must not be under treatment by another endocrinologist. Your Onduo physicians and pharmacists may communicate with your healthcare provider to improve your quality of care.

17.6 Diabetes device data dashboard

As part of the Onduo Type 1 Diabetes program, you'll have access to a third party data dashboard. Your data from any insulin pumps, insulin pens, continuous glucose monitors and blood glucose

meters that you use will be collected there. You can review this data at any time, and can download and print reports for clinicians. Your Onduo Care Specialist will help you set up your dashboard, connect your devices, and show you how to use them. Your device information can help you and your care team to review your diabetes status and make lifestyle and medical changes to improve your health.

18.0 Settings

Customize your Onduo experience through the Settings menu and adjust your schedule, medications, time in range, target range, step goal, emergency contact information, and notifications. You can also connect your meter or check the status of your meter in the app Settings.

18.1 Set and change your schedule

To get meal logging and medication logging notifications at the right time, you can set and change your schedule.

1. Open the Onduo app.
2. Tap **Settings**  > **Schedule** .
3. Tap a part of your schedule to change the time.
4. Tap **Save** when done.

18.2 Add and confirm medications

If your Care Team sets a medication logging goal, you can get helpful notifications on when to log your medications by adding them to the app.

Refer to **Section 9.1** for how to add and confirm medications.

18.3 Adjust your step goal

Track your steps in the Onduo app to better understand how your movement impacts your glucose readings. This information is visible to both you and your Care Team—and you can work together to set goals for yourself.

Refer to **Section 10.1** for how to track steps and adjust your step goal.

18.4 (Old flow - No longer in use) My Appointments

This flow will no longer be used going forward, but instructions are provided for reference for members who had it triggered previously and have not yet completed it.

You can view your upcoming scheduled appointments.

1. Open the Onduo app.
2. Tap **Settings**  > **My appointments**.
3. You will be able to see the following appointments:
 - Office visits with your Health Care Provider that you have told us about
 - Video Visits with Onduo Care Team (e.g.: Onduo MD, licensed mental health care providers (e.g.: psychiatrist or therapist). Cancel or reschedule these appointments using the pencil icon.

18.5 Add and view connected meters

You can add and view your connected meters from the Settings menu.

1. Open the Onduo app.
2. Tap **Settings**  > **My Devices**  .
3. To connect to a disconnected device, tap on the device and click "Pair device".
4. To add a new device, click the "Add device button" at the bottom of the "My devices" page.
5. Devices that are connected will show under the "Connected devices" section, and the disconnected ones under "Disconnected devices".

18.6 Set up and manage notifications

You are automatically reminded when to log information in the Onduo app. All notifications are on by default. However, you can adjust the type of notifications you receive.

1. Tap **Settings**  > **Notifications**  .
2. Tap to toggle a notification type on or off.
3. When done, tap Save or Close  to discard your changes and return to the Home tab.

If you are an iPhone user, be sure to set your device lock and device notification privacy settings appropriately.

18.7 Change your target range and time in range (Type 2 diabetes only)

Important: Always consult with a physician before changing your target range.

If you are assigned a Care Lead, your Care Lead will see the change to your target range or time in range.

18.7.1 Change your target range (Type 2 diabetes only)

The American Diabetes Association (ADA) recommends that you maintain a glucose range of 80–180 mg/dL. You can adjust your target range.

Refer to **Section 6.3.1** for how to adjust your target range.

18.7.2 Change your time in range (Type 2 diabetes only)

If you use a CGM and it seems either too difficult or too easy to get any readings in your target range for the set amount of time, you can manually adjust your time in range goal. If you are assigned a Care Lead, your Care Lead will see the change to your time in range goal.

Refer to **Section 6.3.1** for how to adjust your time in range.

18.8 Change the time (Type 2 diabetes only)

It's important to know when you logged your glucose readings over the course of the day. By default, the app uses your phone's time and date when logging a glucose reading.

If your app displays the wrong time, you can try changing your phone's date and time.

18.9 Add a profile photo

To help your Care Team get to know you, upload or add a profile photo of yourself or something that you feel represents you:

1. Open the Onduo app.
2. Tap **Settings**  > **Edit**.
3. From the menu, choose how you'd like to add a photo:
 - **Take new photo:** Take a photo of yourself using your phone's camera.
 - **Upload photo:** Add a photo of yourself from your phone's photo library (recommended).
 - **Choose avatar from library:** Pick an avatar that represents you. Tap Close  in the top left when done.

18.10 Edit emergency contact information


You can update your emergency contact information in Settings. It's important to keep this up to date, since your Care Team may contact them in case of an emergency. (Note, however, that Onduo is not intended to provide urgent response in the case of an emergency, so if you are having an emergency, you should contact 911 or go to an emergency room.)

1. Open the Onduo app.
2. Tap **Settings** .
3. Tap  **Emergency contacts**
4. Make any necessary edits, then click **Save**. To cancel, click **Cancel**.


18.11 Access your Care Team contact information (for OSA members only)

You can find contact information for your obstructive sleep apnea physician group and Durable Medical Equipment (DME) supplier.

1. Open to Onduo App.
2. Tap **Settings** .
3. Tap **My care team**.

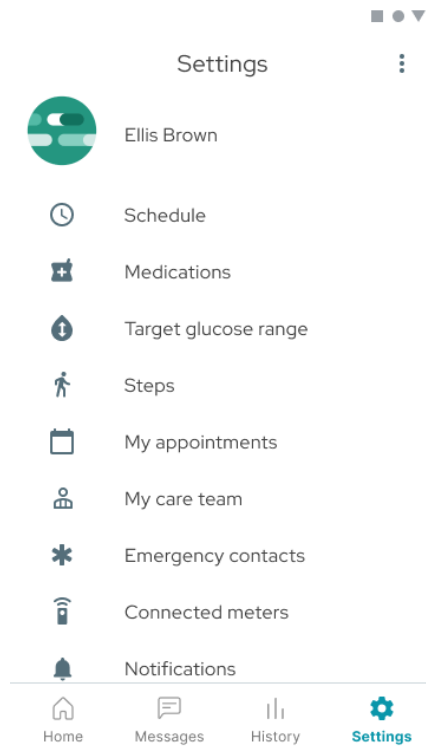Settings screen

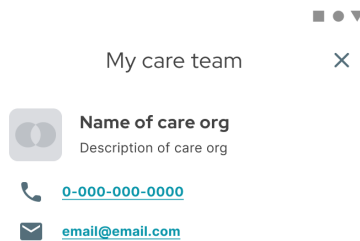Sample Care team view

18.12 Generate clinician report

You can generate a clinician report to take with you on your next doctor's appointment. This report will display a consolidated view of your historical data in PDF format that you can then save, send electronically, download to your device, or print. Please note that if you're a Type 1 diabetes member, you can receive guidance from your T1D console member on generating a clinician report.

1. Open the Onduo app.
2. Tap **Settings** .
3. Tap **For my doctor - BGM** or **For my doctor - CGM** or **For my doctor - Blood Pressure Clinician Report**.

18.13 Update Privacy Settings

18.13.1 Location sharing

You can update your preference to share your location in the background to show health insights based on your location and activity.

1. Open the Onduo app.
2. Tap **Settings** .
3. Tap **Privacy**.
4. Tap to toggle sharing your location on or off.

18.13.2 Import prescription data

You can update your preference for importing your prescription data into the app. If you toggle to off, the app will not import any new medications. Existing medications in your list will not be removed.

1. Open the Onduo app.
2. Tap **Settings** .
3. Tap **Privacy**.
4. Tap to toggle prescription imports on or off.

18.14 Sign out of the app

When you sign out of the app, you will no longer receive notifications or be able to message your Care Team.

1. Open the Onduo app.
2. Tap **Settings** .
3. Tap the Menu  > **Sign out**.

Signing out of the app does not remove your data from the Onduo app.

PRECAUTION

If you have not enabled a passcode on your phone, you will be automatically logged out of the Onduo app after ten minutes of inactivity in order to keep your data confidential. If you are logged out, you'll be directed to your phone's settings to set up a passcode.

18.15 Change password

If you created an Onduo account using your email during registration, you can change the password associated with your account.

1. Open the Onduo app.
2. Tap **Settings**  > **Change password**.
3. Follow the on screen instructions including:
 - a. Entering your current password
 - b. Choosing a new password

19.0 Send feedback about the app

Tell us about your experience using the Onduo app, including what you like, what you don't like, or about issues you've found. Your feedback will be for used product improvement and will not be seen by your Care Team. You may not receive a reply.

You can also call or email support at support@onduo.com or 1-833-HiOnduo (1-833-446-6386). If you're participating in a clinical study, please contact the study support staff.

Send feedback:

1. Open the Onduo app.
2. Tap **Settings** .
3. In the top right, tap Menu  > **Send feedback**.
4. An email will appear with some pre-populated information.
If you receive an error, try copying the email address into your primary email app.
5. Fill out the email with your feedback, then send. Please do not include health information in this feedback.

You should receive a confirmation email to the email account you used to send this feedback. If you do not, call 1-833-HiOnduo (1-833-446-6386) to report your feedback directly. If you're participating in a clinical study, please contact the study support staff.

20.0 Troubleshooting

Learn how to fix common issues with the Onduo app.

20.1 Recover your password

When you signed up for Onduo, if you used a Google Account, you can recover your password using Google's account recovery process.

When you signed up for Onduo, if you created a new account using your email, you can recover your password by clicking on the forgot password link in the sign in with the email screen.

20.2 Fix problems with your BGM, CGM, or Blood Pressure Monitor

If you're having issues with your blood glucose meter (BGM), refer to the meter's instructions for use. If you're having issues with your continuous glucose monitor (CGM), refer to the CGM's instructions for use. If you're having issues with your blood pressure monitor, refer to the blood pressure monitor's instructions for use. If you continue to experience issues syncing your device, please contact Support (see [Section 20.4 Contact support](#)) and inform them of the device brand and model.

20.3 Resolve problems with calling your Care Team

If your Care Team hasn't set their phone number yet, you might receive an error message when trying to call. In this case, message your Care Team and let them know that you'd like to set up a call.

20.4 Contact support

If you're having problems with the app, you can contact the Onduo support team. Call or email support at support@onduo.com or 1-833-HiOnduo (1-833-446-6386). If you're participating in a clinical study, please contact the study support staff.

Additionally, you can contact Onduo support through the app. To contact support:

1. Open the Onduo app.

2. Tap **Settings** .
3. In the top right, tap Menu  > **Contact support**
4. An email will appear with some pre-populated information.
If you receive an error, try copying the email address into your primary email app.
5. Fill out the email with your request for support. You will be asked to grant permission for your device logs to be included with the email. Device logs contain information on how your app and device are working, and what's going wrong. These logs help Onduo developers solve problems faster. Your device logs are encrypted and will never contain your personal information.
6. You should receive a confirmation email to the email account used to send this feedback. If you do not, call 1-833-HiOnduo (1-833-446-6386) to request support directly. If you're participating in a clinical study, please contact the study support staff.

Manufacturer

Verily Life Sciences LLC

999 Bayhill Drive

San Bruno, CA 94066

Trademarks

“VERILY” is a trademark of Verily Life Sciences LLC. All rights reserved.

“Onduo” is a trademark of Onduo LLC. All rights reserved.

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“VERILY” es una marca registrada de Verily Life Sciences LLC. Todos los derechos reservados.
“Onduo” es una marca registrada de Onduo LLC. Todos los derechos reservados.